

Product datasheet

Anti-C1orf114 antibody ab121934

画像数 2

製品の概要

| | |
|--------------|---|
| 製品名 | Anti-C1orf114 antibody |
| 製品の詳細 | Rabbit polyclonal to C1orf114 |
| 由来種 | Rabbit |
| アプリケーション | 適用あり: ICC/IF, IHC-P |
| 種交差性 | 交差種: Human |
| 免疫原 | antigen sequence: KSDASIIEMA CEKEENINQD LKENETVMEH TKRHSDPKS LQDEVSPRRN DIISVPGIQP LDPISDSDSE NSFQESKL, corresponding to amino acids 32-109 of Human C1orf114. Run BLAST with ExPASy Run BLAST with NCBI |
| ポジティブ・コントロール | Human fallopian tube tissue. |

製品の特性

| | |
|--------|---|
| 製品の状態 | Liquid |
| 保存方法 | Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles. |
| バッファー | pH: 7.20 Preservative: 0.02% Sodium azide Constituents: 59% PBS, 40% Glycerol |
| 精製度 | Immunogen affinity purified |
| ポリ/モノ | ポリクローナル |
| アイソタイプ | IgG |

アプリケーション

Our [Abpromise guarantee](#) covers the use of **ab121934** in the following tested applications.


The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

| アプリケーション | Abreviews | 特記事項 |
|----------|-----------|------|
| | | |

| アプリケーション | Abreviews | 特記事項 |
|----------|-----------|---|
| ICC/IF | | Use a concentration of 1 - 4 µg/ml. Recommend PFA Fixation and Triton X-100 treatment |
| IHC-P | | 1/900. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol. |


ターゲット情報

画像

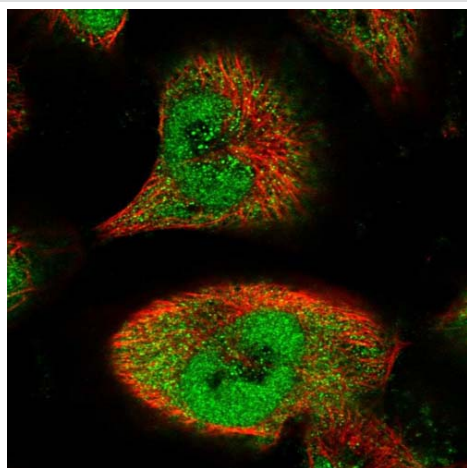

Immunofluorescent staining of Human cell line U-251 MG shows positivity in nucleus but not nucleoli and cytoplasm. Recommended concentration of ab121934 1-4 µg/ml. Cells treated with PFA/Triton X-100.

Immunocytochemistry/ Immunofluorescence - Anti-C1orf114 antibody (ab121934)

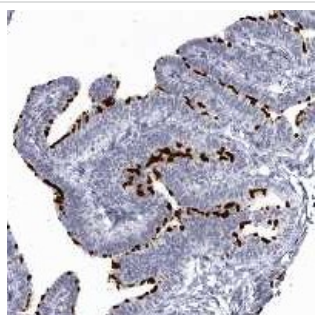

ab121934, at 1/900 dilution, staining C1orf114 in paraffin embedded Human fallopian tube tissue by Immunohistochemistry.

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-C1orf114 antibody (ab121934)

Please note: All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours

- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <http://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors