

Product datasheet

Recombinant human PI 3 Kinase p110 alpha + PI 3 kinase p85 alpha protein ab90730

画像数 2

カルタヘナ法

製品の概要

製品名 Recombinant human PI 3 Kinase p110 alpha + PI 3 kinase p85 alpha protein
タンパク質長 Full length protein

製品の詳細

由来 Recombinant
由来 Baculovirus infected Sf9 cells
アミノ酸配列
生物種 Human

特性

Our [Abpromise guarantee](#) covers the use of **ab90730** in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

生理活性 Activity: 9,900 units/mg 1 unit equals 1 picomole phosphate transferred to PIP2 per min. Activity test: ATP depletion assay, PI as substrate, 1 μM ATP.

アプリケーション Functional Studies
SDS-PAGE

精製度 > 90 % SDS-PAGE.

製品の状態 Liquid

前処理および保存


保存方法および安定性 Shipped on dry ice. Upon delivery aliquot and store at -80°C. Avoid freeze / thaw cycles.
pH: 7.90
Constituents: 0.05% Brij, 0.79% Tris HCl, 10% Glycerol, 1.74% Sodium chloride
This product is an active protein and may elicit a biological response in vivo, handle with caution.

関連情報

関連性


Phosphoinositide 3-kinases (PI3Ks) phosphorylate phosphatidylinositols (PIs) at their 3' OH position generating lipid second messengers and thereby regulate numerous biological processes including cell growth, differentiation, survival, proliferation, migration and metabolism. On the basis of structural similarities and substrate specificity, the PI3K family can be subdivided into three classes termed I, II, and III. All human class I members are heterodimers consisting of a catalytic subunit (MW approx. 110 kDa) and a non-catalytic subunit (MW 50, 55, 85, or 101 kDa) and are known to phosphorylate phosphatidylinositol (PI), phosphatidylinositol-4-mono-phosphate (PIP) and phosphatidylinositol-4,5-bisphosphate (PIP₂) in vitro. The class I members can be further subdivided into class IA and IB PI3Ks. Class IA exists in three isoforms (p110alpha, p110beta and p110delta whereas the only class IB member is termed p110gamma. Class IA PI3Ks are activated by adaptor proteins such as Ras or BCAP, or tyrosine-kinase-associated receptors including antigen, co-stimulatory and cytokine receptors (e.g. CD19, CD28, Insulin receptor, EGFR, and PDGFR). p110gamma is activated by G-protein-coupled receptors (GPCRs). Effectors of class I PI3Ks are pleckstrin homology domain proteins such as Akt/PKB, BTK, TEC, ITK, BAM32, and small GTPases (e.g. Cdc42, Rac, or Ras). The action of PI3Ks is regulated by the phosphatidylinositol-3,4,5-trisphosphate phosphatases SHIP and PTEN.

画像

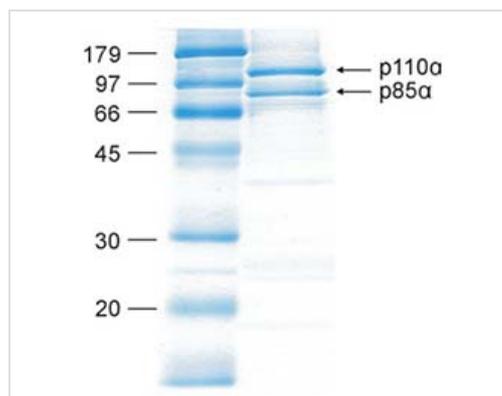

SDS-PAGE analysis of 3µg of ab90730

SDS-PAGE - Recombinant human PI 3 Kinase p110 alpha + PI 3 kinase p85 alpha protein (ab90730)

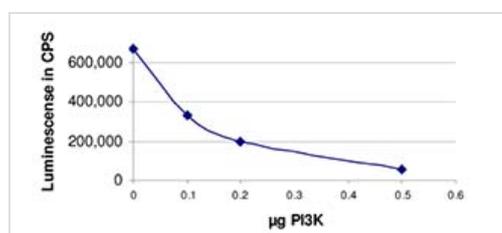

ATP depletion assay, PI as substrate, 1 µM ATP

Functional Studies - Recombinant human PI 3 Kinase p110 alpha + PI 3 kinase p85 alpha protein (ab90730)

Please note: All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours

- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors