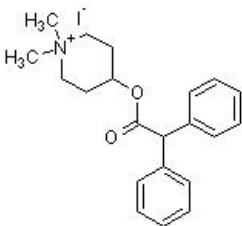


4-DAMP, Muscarinic receptor antagonist ab120144


4 References 画像数 2

製品の概要


製品名	4-DAMP, Muscarinic receptor antagonist
製品の詳細	Muscarinic receptor antagonist
生理活性の詳細	Muscarinic antagonist. (K_i values for human receptors are 0.57 (M_1), 7.3 (M_2), 0.37 (M_3), 0.72 (M_4) and 0.55 nM (M_5)).
CAS 番号	1952-15-4
構造式	

製品の特性

体系名	4-Diphenylacetoxy- <i>N</i> -methylpiperidine methiodide
分子量	451.35
分子式	$C_{21}H_{26}INO_2$
PubChem 登録番号	3014059
保存方法	Store at Room Temperature. The product can be stored for up to 12 months.
溶解性	Soluble in DMSO to 100 mM and in ethanol to 25 mM
使用に関する注意	<p>Wherever possible, you should prepare and use solutions on the same day. However, if you need to make up stock solutions in advance, we recommend that you store the solution as aliquots in tightly sealed vials at -20°C. Generally, these will be useable for up to one month. Before use, and prior to opening the vial we recommend that you allow your product to equilibrate to room temperature for at least 1 hour.</p> <p>Need more advice on solubility, usage and handling? Please visit our frequently asked questions (FAQ) page for more details.</p>
SMILES 線形表記	<chem>[I-].O=C(OC1CC[N+](C)(C)CC1)C(c2ccccc2)c3ccccc3</chem>
由来	Synthetic

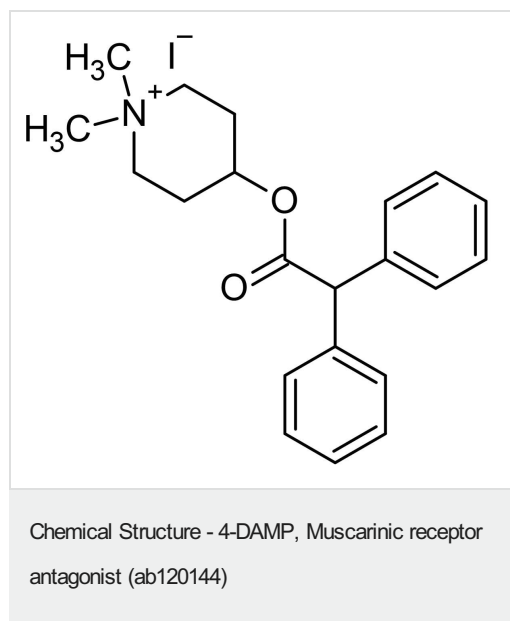

2D chemical structure image of ab120144, 4-DAMP, Muscarinic receptor antagonist

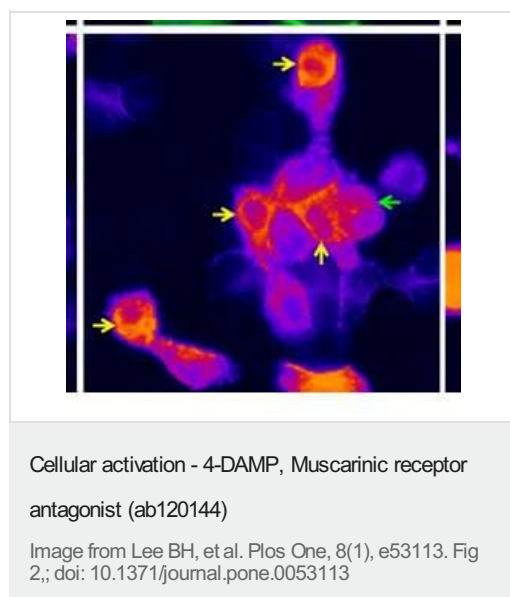

(A) rhAQP5 vector-transfected HSG cells were stimulated with CCh (100 μ M) in the presence or absence of 4-DAMP (10 μ M), prenzepine (10 μ M), and mixture of 4-DAMP and prenzepine (M3R or M1R specific antagonist, respectively).

(B) To monitor hAQP5 trafficking, original GFP pictures in (A) were converted into 'Fire' 3D LUT (lookup tables) setting in ImageJ software and only GFP signal above the fluorescent signal threshold (?80, red or yellow color in converted picture), was captured in (B) as positive GFP signal. Yellow arrows indicate the positive GFP signal on the membrane and green arrows indicates signal in the cytoplasm in HSG cells.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES, NOT FOR USE IN HUMANS"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors
- Abcam biochemicals are novel compounds and we have not tested their biological activity in house. Please use the literature to identify how to use these products effectively. If you require further assistance please contact the scientific support team