

### Anti-ACRBP antibody ab64809

★★★★☆ [1 Abreviews](#) [4 References](#) [画像数 2](#)

#### 製品の概要

製品名	Anti-ACRBP antibody
製品の詳細	Rabbit polyclonal to ACRBP
由来種	Rabbit
アプリケーション	<b>適用あり:</b> WB, IHC-P
種交差性	<b>交差種:</b> Human
免疫原	A synthetic peptide derived from an internal region of human ACRBP.
特記事項	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&amp;As</p>

#### 製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at -20°C. Stable for 12 months at -20°C.
バッファー	<p>pH: 7.40</p> <p>Preservative: 0.02% Sodium azide</p> <p>Constituents: PBS, 50% Glycerol, 0.87% Sodium chloride</p>
精製度	Immunogen affinity purified
ポリ/モノ	ポリクローナル
アイソタイプ	IgG

#### アプリケーション


**The Abpromise guarantee**      **Abpromise保証は、次のテスト済みアプリケーションにおけるab64809の使用に適用されます**  
 アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
WB	★★★★★ (1)	1/500 - 1/1000. Detects a band of approximately 61 kDa (predicted molecular weight: 61 kDa).
IHC-P		Use a concentration of 1 µg/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

## ターゲット情報


機能	May be involved in packaging and condensation of the acrosin zymogen in the acrosomal matrix via its association with proacrosin.
組織特異性	Expression restricted to testis in normal tissue. Expressed in a wide spectrum of cancers, including bladder, breast, liver, lung and colon cancers.
翻訳後修飾	The N-terminus is blocked. Phosphorylated on Tyr residues in capacitated sperm.
細胞内局在	Secreted. Cytoplasmic vesicle, secretory vesicle, acrosome. Colocalizes with proacrosin in the acrosome of sperm.

## 画像

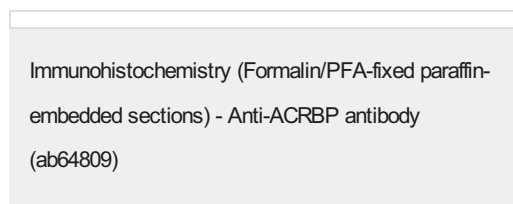

IHC image of ab64809 staining in human breast carcinoma formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab64809, 1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

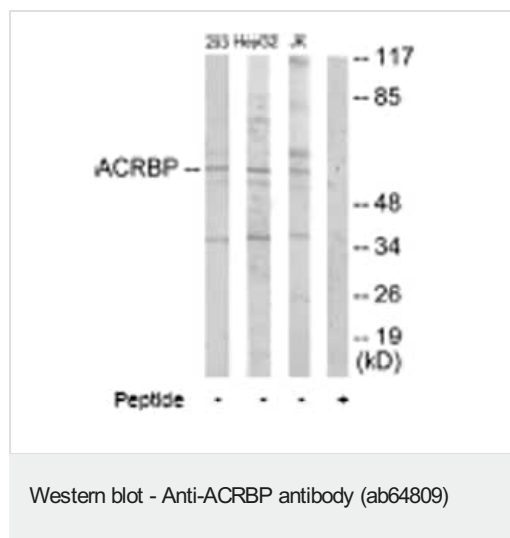

**All lanes :** Anti-ACRBP antibody (ab64809) at 1/500 dilution

**Lane 1 :** Extracts of 293 cells

**Lane 2 :** Extracts of HepG2 cells

**Lane 3 :** Extracts of Jurkat cells

**Lane 4 :** Extracts of Jurkat cells with immunizing Peptide

Lysates/proteins at 5 µg per lane.

**Predicted band size:** 61 kDa

**Observed band size:** 61 kDa

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

### Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

### Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors