abcam

Product datasheet

Alexa Fluor® 488 Anti-ACTN3 antibody [EP2531Y] ab210420

יעלאעבע RabMAb

画像数3

製品の概要

特記事項

製品名 Alexa Fluor® 488 Anti-ACTN3 antibody [EP2531Y]

製品の詳細 Alexa Fluor® 488 Rabbit monoclonal [EP2531Y] to ACTN3

由来種 Rabbit

標識 Alexa Fluor® 488. Ex: 495nm, Em: 519nm

アプリケーション 適用あり: IHC-Fr

種交差性 交差種: Rat

交差が予測される動物種: Mouse, Human 4

免疫原 Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール IHC-Fr: normal rat skeletal muscle tissue sections

> Our RabMAb® technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb**® **patents**.

Alexa Fluor[®] is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor® dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor® dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor® dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or

outlicensing@thermofisher.com.

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

バッファー pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: PBS, 30% Glycerol (glycerin, glycerine), 1% BSA

精製度 Affinity purified

ポリ/モノ モノクローナル

クローン名 EP2531Y

アイソタイプ lgG

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab210420の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-Fr		1/100.

ターゲット情報

機能 F-actin cross-linking protein which is thought to anchor actin to a variety of intracellular structures.

This is a bundling protein.

組織特異性 Expressed only in a subset of type 2 skeletal muscle fibers.


配列類似性 Belongs to the alpha-actinin family.

Contains 1 actin-binding domain.


Contains 2 CH (calponin-homology) domains.

Contains 2 EF-hand domains. Contains 4 spectrin repeats.

画像

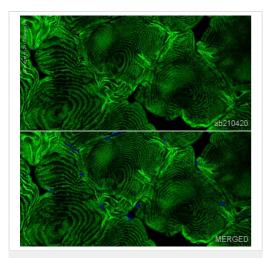

Immunohistochemistry (Frozen sections) - Alexa Fluor® 488 Anti-ACTN3 antibody [EP2531Y] (ab210420)

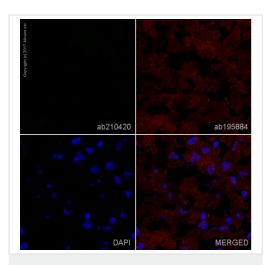

Immunohistochemistry (Frozen sections) - Alexa Fluor® 488 Anti-ACTN3 antibody [EP2531Y] (ab210420)

IHC image of ACTN3 staining in a section of frozen normal rat skeletal muscle.

The section was fixed using 10% formaldehyde in 1XPBS for 10 minutes. No antigen retrieval step was performed prior to staining. Non-specific protein-protein interactions were then blocked in TBS containing 0.025% (v/v) Triton X-100, 0.3M (w/v) glycine and 1% (w/v) BSA for 1h at room temperature. The section was then incubated overnight at +4°C in TBS containing 0.025% (v/v) Triton X-100 and 1% (w/v) BSA with ab210420 at 1/100 dilution (shown in green). Nuclear DNA was labelled with DAPI (shown in blue). The section was then mounted using Fluoromount[®].

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).


For other IHC staining systems (automated and non-automated), customers should optimize variable parameters such as antigen retrieval conditions, antibody concentrations and incubation times.

Negative IHC image of ACTN3 staining in a section of frozen normal rat liver.

The section was fixed using 10% formaldehyde in 1XPBS for 10 minutes. No antigen retrieval step was performed prior to staining. Non-specific protein-protein interactions were then blocked in TBS containing 0.025% (v/v) Triton X-100, 0.3M (w/v) glycine and 1% (w/v) BSA for 1h at room temperature. The section was then incubated overnight at +4°C in TBS containing 0.025% (v/v) Triton X-100 and 1% (w/v) BSA with ab210420 at 1/100 dilution (shown in green) and counterstained using ab195884, Rat monoclonal to Tubulin (Alexa Fluor® 647), at 1/250 dilution (shown in red). Nuclear DNA was labelled with DAPI (shown in blue). The section was then mounted using Fluoromount®.

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

For other IHC staining systems (automated and non-automated), customers should optimize variable parameters such as antigen retrieval conditions, antibody concentrations and incubation times.

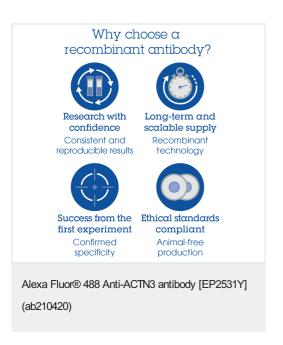

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors