abcam

Product datasheet

Alexa Fluor® 488 Anti-CD34 antibody [EP373Y] ab195013

יעלאעבע RabMAb

★★★★★ 5 Abreviews 4 References 画像数 2

製品の概要

免疫原

製品名 Alexa Fluor® 488 Anti-CD34 antibody [EP373Y]

製品の詳細 Alexa Fluor® 488 Rabbit monoclonal [EP373Y] to CD34

由来種 Rahhit

標識 Alexa Fluor® 488, Ex: 495nm, Em: 519nm

アプリケーション 適用あり: Flow Cyt (Intra)

種交差性 交差種: Human

交差が予測される動物種: Mouse, Sheep, Dog, Pig, African bush elephant 4 非交差種: Rat

Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール Flow Cyt (intra): HUVEC cells.

Our RabMAb® technology is a patented hybridoma-based technology for making rabbit 特記事項 monoclonal antibodies. For details on our patents, please refer to RabMAb® patents.

> Alexa Fluor[®] is a registered trademark of Molecular Probes, Inc., a Thermo Fisher Scientific Company. The Alexa Fluor® dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor® dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor[®] dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or

outlicensing@thermofisher.com.

製品の特性


製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

解離定数(Kn值)

 $K_D = 1.15 \times 10^{-10} M$

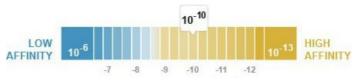

Learn more about K_D

バッファー pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: PBS, 30% Glycerol (glycerin, glycerine), 1% BSA

精製度 Protein A purified

ポリ/モノ モノクローナル

クローン名 EP373Y

アイソタイプ lgG

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab195013の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
Flow Cyt (Intra)		1/500.

ターゲット情報

機能 Possible adhesion molecule with a role in early hematopoiesis by mediating the attachment of

stem cells to the bone marrow extracellular matrix or directly to stromal cells. Could act as a scaffold for the attachment of lineage specific glycans, allowing stem cells to bind to lectins expressed by stromal cells or other marrow components. Presents carbohydrate ligands to

selectins.

組織特異性 Selectively expressed on hematopoietic progenitor cells and the small vessel endothelium of a

variety of tissues.

配列類似性 Belongs to the CD34 family.


発生段階 On early hematopoietic progenitor cells.

翻訳後修飾 Highly glycosylated.

Phosphorylated on serine residues by PKC.

細胞内局在 Membrane.

画像

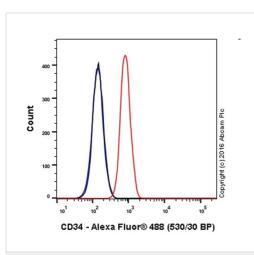


Flow Cytometry (Intracellular) - Alexa Fluor® 488 Anti-CD34 antibody [EP373Y] (ab195013) Overlay histogram showing HUVEC cells stained with ab195013 (red line). The cells were fixed with 4% formaldehyde (10 min) and then permeabilized with 0.1% PBS-Triton X-100 for 15 min. The cells were then incubated in 1x PBS / 10% normal goat serum to block non-specific protein-protein interactions followed by the antibody (ab195013, 1/500 dilution) for 30 min at 22°C.

Isotype control antibody (black line) was Rabbit IgG (monoclonal) Alexa Fluor® 488 (<u>ab199091</u>) used at the same concentration and conditions as the primary antibody. Unlabelled sample (blue line) was also used as a control.

Acquisition of >5,000 events were collected using a 50 mW Blue laser (488nm) and 530/30 bandpass filter.

This antibody gave a positive signal in HUVEC cells fixed with 80% methanol (5 min)/permeabilized with 0.1% PBS-Triton X-100 for 15 min used under the same conditions.

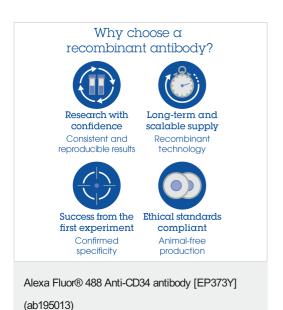

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you

• We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors