

# Anti-Alpha 1 microglobulin antibody ab47980

画像数 2

### 製品の概要

製品名	Anti-Alpha 1 microglobulin antibody
製品の詳細	Rabbit polyclonal to Alpha 1 microglobulin
由来種	Rabbit
アプリケーション	<b>適用あり:</b> ICC/IF, IHC-P
種交差性	<b>交差種:</b> Human
免疫原	Human Alpha 1 Microglobulin
特記事項	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&amp;As</p>

### 製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long term.
バッファー	<p>pH: 7.40</p> <p>Preservative: 0.02% Sodium azide</p> <p>Constituents: 49.98% PBS, 50% Glycerol (glycerin, glycerine)</p>
精製度	Protein G purified
ポリ/モノ	ポリクローナル
アイソタイプ	IgG


### アプリケーション

**The Abpromise guarantee**      **Abpromise保証は、** 次のテスト済みアプリケーションにおけるab47980の使用に適用されます  
 アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
ICC/IF		Use a concentration of 1 - 5 µg/ml.
IHC-P		Use a concentration of 4 µg/ml.


ターゲット情報	
機能	Inter-alpha-trypsin inhibitor inhibits trypsin, plasmin, and lysosomal granulocytic elastase. Inhibits calcium oxalate crystallization. Trypstatin is a trypsin inhibitor.
組織特異性	Expressed by the liver and secreted in plasma. Alpha-1-microglobulin occurs in many physiological fluids including plasma, urine, and cerebrospinal fluid. Inter-alpha-trypsin inhibitor is present in plasma and urine.
配列類似性	In the N-terminal section; belongs to the calycin superfamily. Lipocalin family. Contains 2 BPTI/Kunitz inhibitor domains.
翻訳後修飾	The precursor is proteolytically processed into separately functioning proteins. 3-hydroxykynurenine, an oxidized tryptophan metabolite that is common in biological fluids, reacts with Cys-53, Lys-111, Lys-137, and Lys-149 to form heterogeneous polycyclic chromophores including hydroxanthommatin. The reaction by alpha-1-microglobulin is autocatalytic; the human protein forms chromophore even when expressed in insect and bacterial cells. The chromophore can react with accessible cysteines forming non-reducible thioether cross-links with other molecules of alpha-1-microglobulin or with other proteins such as Ig alpha-1 chain C region 'Cys-352'. Heavy chains are interlinked with bikunin via a chondroitin 4-sulfate bridge to the their C-terminal aspartate. Addition of glycosaminoglycan chondroitin sulfate, allows cross-linking between the different components.
細胞内局在	Secreted.

画像

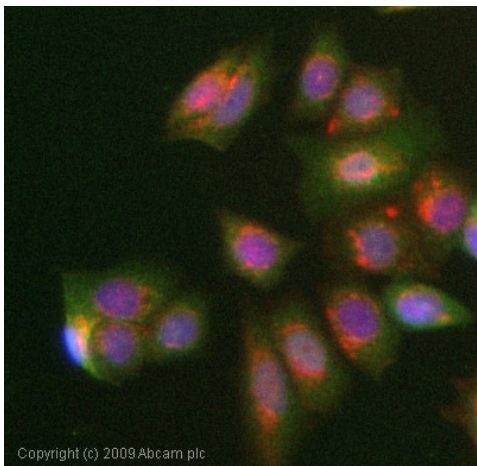

Immunocytochemistry/ Immunofluorescence - Anti-Alpha 1 microglobulin antibody (ab47980)

ICC/IF image of ab47980 stained MCF7 cells. The cells were 4% PFA fixed (10 min) and then incubated in 1%BSA / 10% normal goat serum / 0.3M glycine in 0.1% PBS-Tween for 1h to permeabilise the cells and block non-specific protein-protein interactions. The cells were then incubated with the antibody (ab47980, 1µg/ml) overnight at +4°C. The secondary antibody (green) was Alexa Fluor® 488 goat anti-rabbit IgG (H+L) used at a 1/1000 dilution for 1h. Alexa Fluor® 594 WGA was used to label plasma membranes (red) at a 1/200 dilution for 1h. DAPI was used to stain the cell nuclei (blue) at a concentration of 1.43µM.

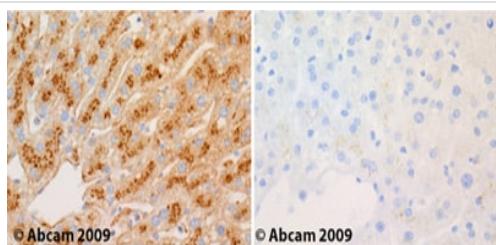

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Alpha 1 microglobulin antibody (ab47980)

ab47980 Microglobulin alpha 1 in human liver.

Left panel: with primary antibody at 4 µg/ml. Right panel: isotype control.

Sections were stained using an automated system (DAKO Autostainer Plus), at room temperature: sections were rehydrated and antigen retrieved with the Dako 3 in 1 AR buffers EDTA pH 9.0. Slides were peroxidase blocked in 3% H<sub>2</sub>O<sub>2</sub> in methanol for 10 mins. They were then blocked with Dako Protein block for 10 minutes (containing casein 0.25% in PBS) then incubated with primary antibody for 20 min and detected with Dako envision flex amplification kit for 30 minutes. Colorimetric detection was completed with Diaminobenzidine for 5 minutes. Slides were counterstained with Haematoxylin and coverslipped under DePeX. Please note that for manual staining we recommend to optimize the primary antibody concentration and incubation time (overnight incubation), and amplification may be required.

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

### Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise,

please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

## **Terms and conditions**

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors