abcam

Product datasheet

Anti-Apelin antibody ab59469

18 References 画像数 1

製品の概要

免疫原

製品名 Anti-Apelin antibody

製品の詳細 Rabbit polyclonal to Apelin

由来種 Rabbit

特異性 This antibody reacts specifically with 8,5 kDa of human APLN protein.

アプリケーション 適用あり: IHC-P

種交差性 交差種: Human

交差が予測される動物種: Mouse, Rat, Cow 🕰

Synthetic peptide corresponding to Human Apelin (C terminal).

Database link: Q9ULZ1

特記事項The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long

term.

バッファー Constituent: Whole serum

精製度 Whole antiserum ポリ/モノ ポリクローナル

アイソタイプ lgG

アプリケーション


アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		1/500.

ターゲット情報

機能	Endogenous ligand for APJ, an alternative coreceptor with CD4 for HIV-1 infection. Inhibits HIV-1 entry in cells coexpressing CD4 and APJ. Apelin-36 has a greater inhibitory activity on HIV infection than other synthetic apelin derivatives. The oral intake in the colostrum and the milk could have a role in the modulation of the immune responses in neonates. May also have a role in the central control of body fluid homeostasis by influencing AVP release and drinking behavior.
組織特異性	Expressed in the brain with highest levels in the frontal cortex, thalamus, hypothalamus and midbrain. Secreted by the mammary gland into the colostrum and the milk.
配列類似性	Belongs to the apelin family.
翻訳後修飾	Several active peptides may be produced by proteolytic processing of the peptide precursor.
細胞内局在	Secreted.

画像

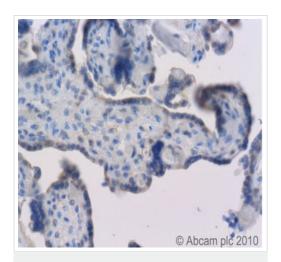

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Apelin antibody (ab59469)

ab59469 (1:500) staining Apelin in human placenta using an automated system (DAKO Autostainer Plus). Using this protocol there is strong staining of cytoplasmic/secreted regions of cells in the basal uterine glands.

Sections were rehydrated and antigen retrieved with the Dako 3 in 1 AR buffer EDTA pH 9.0 in a DAKO PT link. Slides were peroxidase blocked in 3% H2O2 in methanol for 10 mins. They were then blocked with Dako Protein block for 10 minutes (containing casein 0.25% in PBS) then incubated with primary antibody for 20 min and detected with Dako envision flex amplification kit for 30 minutes. Colorimetric detection was completed with Diaminobenzidine for 5 minutes. Slides were counterstained with Haematoxylin and coverslipped under DePeX. Please note that, for manual staining, optimization of primary antibody concentration and incubation time is recommended. Signal amplification may be required.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

• Replacement or refund for products not performing as stated on the datasheet

- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors