abcam

Product datasheet

Anti-CEACAM5 + CEACAM6 antibody [NCRC16 (AKA 161)] ab22705

2 Abreviews 1 References 画像数 2

製品の概要

製品名 Anti-CEACAM5 + CEACAM6 antibody [NCRC16 (AKA 161)]

製品の詳細 Mouse monoclonal [NCRC16 (AKA 161)] to CEACAM5 + CEACAM6

由来種 Mouse

特異性 Cross reacts with a similar epitope on NCA.

アプリケーション 適用あり: ICC/IF, IHC-P

種交差性 交差種: Human

免疫原 Colorectal cancer cell lines (Human).

ポジティブ・コントロール LS174 T cells, colon carcinoma IHC-P: FFPE human colon tissue sections.

特記事項

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

רק"ע". Preservative: 0.02% Sodium azide

Constituent: 99.98% PBS

精製度 Protein A/G purified

ポリ/モノ モノクローナル

クローン名 NCRC16 (AKA 161)

アイソタイプ lgG1

1

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab22705の使用に適用されますアプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。


アプリケーション	Abreviews	特記事項
ICC/IF		Use a concentration of 1 µg/ml.
IHC-P		Use a concentration of 0.1 μ g/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

ターゲット情報

細胞内局在

CEACAM5: Cell membrane. CEACAM6: Cell membrane.


画像

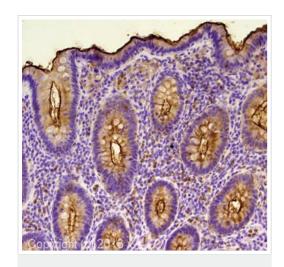

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Carcino Embryonic
Antigen CEA antibody [NCRC16 (AKA 161)]
(ab22705)

IHC image of Carcino Embryonic Antigen CEA staining in human colon formalin fixed paraffin embedded tissue section, performed on a Leica Bond system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab22705, 0.5µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

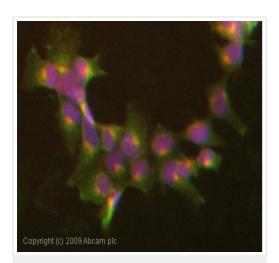

Immunocytochemistry/ Immunofluorescence - Anti-Carcino Embryonic Antigen CEA antibody [NCRC16 (AKA 161)] (ab22705)

ICC/IF image of ab22705 stained Hek293 cells. The cells were 4% formaldehyde fixed (10 min) and then incubated in 1%BSA / 10% normal goat serum / 0.3M glycine in 0.1% PBS-Tween for 1h to permeabilise the cells and block non-specific protein-protein interactions. The cells were then incubated with the antibody (ab22705, 1µg/ml) overnight at +4°C. The secondary antibody (green) was Alexa Fluor® 488 goat anti-mouse IgG (H+L) used at a 1/1000 dilution for 1h. Alexa Fluor® 594 WGA was used to label plasma membranes (red) at a 1/200 dilution for 1h. DAPI was used to stain the cell nuclei (blue) at a concentration of 1.43µM.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors