

Anti-Ganglioside GD3 antibody [R24] ab11779

22 References **画像数 3**

製品の概要

製品名	Anti-Ganglioside GD3 antibody [R24]
製品の詳細	Mouse monoclonal [R24] to Ganglioside GD3
由来種	Mouse
アプリケーション	適用あり: IHC-P, ICC/IF, IHC-Fr
種交差性	交差種: Human
免疫原	Tissue, cells or virus corresponding to Human Ganglioside GD3. This antibody was raised against SK-MEL-28 melanoma cell line.
ポジティブ・コントロール	ICC/IF: SH-SY5Y cells. IHC-P: Human melanoma tissue. IHC-Fr: Human melanoma tissue.
特記事項	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As</p>

製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Store at -20°C or -80°C. Avoid freeze / thaw cycle.
バッファー	Preservative: 0.09% Sodium azide Constituents: 98.9% PBS, 1% BSA
精製度	Protein A purified
特記事項(精製)	Affinity purified from tissue culture supernatant.
ポリ/モノ	モノクローナル
クローン名	R24
アイソタイプ	IgG3
軽鎖の種類	kappa

アプリケーション


The Abpromise guarantee Abpromise保証は、次のテスト済みアプリケーションにおけるab11779の使用に適用されます
アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご確認ください。

アプリケーション	Abreviews	特記事項
IHC-P		Use a concentration of 1 µg/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.
ICC/IF		Use a concentration of 5 µg/ml.
IHC-Fr		1/40.

ターゲット情報

関連性 Almost all melanomas, astrocytomas, a proportion of sarcomas, a small number of carcinomas, some nevi, as well as normal melanocytes express GD3 antigen. It is one of the most important markers for malignant melanoma. Antibodies to the GD3 ganglioside can induce partial remission of tumor growth in animals as well as in Humans via enhancement of cytotoxic and proliferative response of lymphocytes.


画像

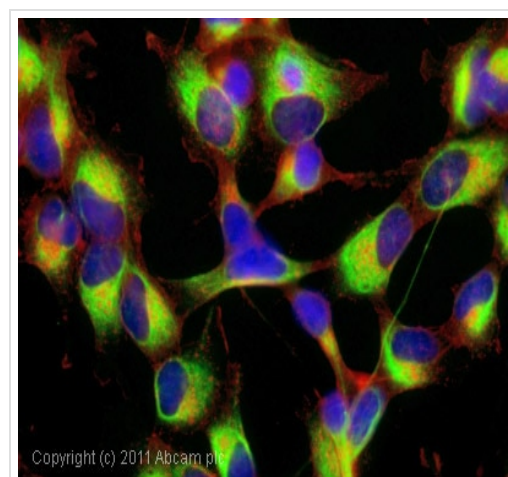

Copyright (c) 2011 Abcam plc


Immunocytochemistry/ Immunofluorescence - Anti-Ganglioside GD3 antibody [R24] (ab11779)

ICC/IF image of ab11779 stained SHSY5Y cells. The cells were 4% formaldehyde fixed (10 min) and then incubated in 1%BSA / 10% normal goat serum / 0.3M glycine in 0.1% PBS-Tween for 1h to permeabilise the cells and block non-specific protein-protein interactions. The cells were then incubated with the antibody (ab11779, 5µg/ml) overnight at +4°C. The secondary antibody (green) was Alexa Fluor® 488 goat anti-mouse IgG (H+L) used at a 1/1000 dilution for 1h. Alexa Fluor® 594 WGA was used to label plasma membranes (red) at a 1/200 dilution for 1h. DAPI was used to stain the cell nuclei (blue) at a concentration of 1.43µM.

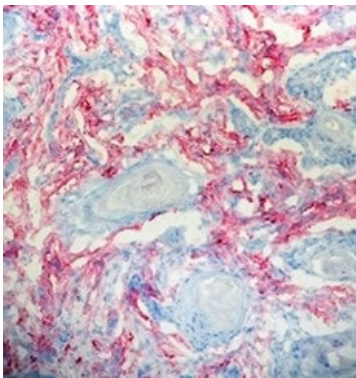

Immunohistochemistry (Frozen sections) - Anti-Ganglioside GD3 antibody [R24] (ab11779)

Frozen section of human melanoma tissue stained for Ganglioside GD3 with ab11779 at 1/40 dilution in immunohistochemical analysis.

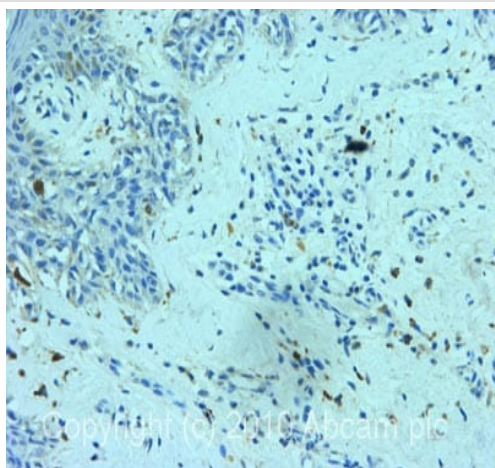

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Ganglioside GD3 antibody [R24] (ab11779)

IHC image of ab11779 staining in human melanoma formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab11779, 1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you

- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors