# abcam

## Product datasheet

## Anti-GDF6 antibody ab73288

★★★★★ 4 Abreviews 2 References 画像数 2

#### 製品の概要

製品名 Anti-GDF6 antibody

製品の詳細 Rabbit polyclonal to GDF6

由来種 Rabbit

アプリケーション **適用あり**: WB

種交差性 交差種: Human

免疫原 A 17 amino acid peptide near the carboxy terminus of human GDF6.

ポジティブ・コントロール SK-N-SH lysate.

特記事項 The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

## 製品の特性

製品の状態 Liquid

**保存方法** Shipped at 4°C. Store at +4°C.

**バッファー** pH: 7.2

Preservative: 0.02% Sodium azide

Constituent: PBS

精製度 Immunogen affinity purified

**ポリ/モノ** ポリクローナル

アイソタイプ lqG

### アプリケーション

**The Abpromise guarantee** <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab73288の使用に適用されますアプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

1

アプリケーション	Abreviews	特記事項
WB	**** <u>(2)</u>	Use a concentration of 0.5 - 1 µg/ml. Detects a band of approximately 55 kDa (predicted molecular weight: 51 kDa).

#### ターゲット情報

#### 機能

Growth factor that controls proliferation and cellular differentiation in the retina and bone formation. Plays a key role in regulating apoptosis during retinal development. Establishes dorsal-ventral positional information in the retina and controls the formation of the retinotectal map (PubMed:23307924). Required for normal formation of bones and joints in the limbs, skull, digits and axial skeleton. Plays a key role in establishing boundaries between skeletal elements during development. Regulation of GDF6 expression seems to be a mechanism for evolving species-specific changes in skeletal strucutres. Seems to positively regulates differentiation of chondrogenic tissue through the growth factor receptors subunits BMPR1A, BMPR1B, BMPR2 and ACVR2A, leading to the activation of SMAD1-SMAD5-SMAD8 complex. The regulation of chondrogenic differentiation is inhibited by NOG (PubMed:26643732). Also involved in the induction of adipogenesis from mesenchymal stem cells. This mechanism acts through the growth factor receptors subunits BMPR1A, BMPR2 and ACVR2A and the activation of SMAD1-SMAD5-SMAD8 complex and MAPK14/p38.

## 関連疾患

Klippel-Feil syndrome 1, autosomal dominant

A chromosomal aberration involving GDF6 has been found in a patient with Klippel-Feil syndrome (KFS). Paracentric inv(8)(q22;2q23.3).

Microphthalmia, isolated, 4

Leber congenital amaurosis 17

Defects in POP1 may be the cause of multiple synostoses syndrome (SYNS). SYNS is a bone disease characterized by multiple progressive joint fusions that commonly involve proximal interphalangeal, tarsal-carpal joints. Additional features can include progressive conductive deafness.


#### 配列類似性

Belongs to the TGF-beta family.

## 細胞内局在

Secreted.

## 画像

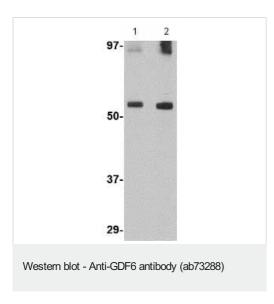


**Lane 1 :** Anti-GDF6 antibody (ab73288) at 0.5 μg/ml **Lane 2 :** Anti-GDF6 antibody (ab73288) at 1 μg/ml

All lanes: SK-N-SH lysate

Lysates/proteins at 15 µg per lane.

**Predicted band size:** 51 kDa **Observed band size:** 55 kDa

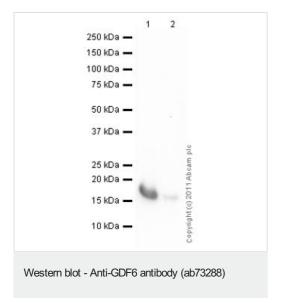

All lanes: Anti-GDF6 antibody (ab73288) at 1 µg/ml

**Lane 1 :** Recombinant Human GDF6 protein (<u>ab50230</u>) at 0.1 μg **Lane 2 :** Recombinant Human GDF6 protein (<u>ab50230</u>) at 0.01 μg

## Secondary

**All lanes :** Goat Anti-Rabbit lgG H&L (HRP) preadsorbed (ab97080) at 1/5000 dilution

Developed using the ECL technique.

Performed under reducing conditions.

Predicted band size: 51 kDa

Exposure time: 30 seconds

<u>ab50230</u> is a homodimer consisting of two 120aa monomers. The homodimer format is expected to run at 27kDa, so the observed band is thought to be the monomer form.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

## Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours

- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

#### Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors