# abcam

# Product datasheet

# Anti-KLB antibody ab106794

11 References 画像数 2

製品の概要

製品名 Anti-KLB antibody

製品の詳細 Rabbit polyclonal to KLB

由来種 Rabbit

特異性 From Mar 2024, QC testing of replenishment batches of this polyclonal changed. All tested and

expected application and reactive species combinations are still covered by our Abcam product promise. However, we no longer test all applications. For more information on a specific batch,

please contact our Scientific Support who will be happy to help.

**アプリケーション 適用あり:** WB, IHC-P

種交差性 交差種: Human

交差が予測される動物種: Macaque monkey, Gorilla 🔒

免疫原 Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール This antibody gave a positive signal in HepG2 whole cell lysate. This antibody gave a positive

result in IHC in the following FFPE tissue: Human normal stomach.

**特記事項**The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

**バッファー** pH: 7.40

Preservative: 0.02% Sodium azide

Constituent: PBS

Batches of this product that have a concentration < 1mg/ml may have BSA added as a stabilising

agent. If you would like information about the formulation of a specific lot, please contact our

scientific support team who will be happy to help.

精製度 Immunogen affinity purified

**ポリ/モノ** ポリクローナル

アイソタイプ lgG

#### アプリケーション

**The Abpromise guarantee** <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab106794の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
WB		Use a concentration of 1 µg/ml. Detects a band of approximately 125 kDa (predicted molecular weight: 120 kDa).
IHC-P		Use a concentration of 5 µg/ml.

# ターゲット情報

機能 Contributes to the transcriptional repression of cholesterol 7-alpha-hydroxylase (CYP7A1), the

rate-limiting enzyme in bile acid synthesis. Probably inactive as a glycosidase. Increases the

ability of FGFR1 and FGFR4 to bind FGF21.

**配列類似性** Belongs to the glycosyl hydrolase 1 family. Klotho subfamily.


ドメイン Contains 2 glycosyl hydrolase 1 regions. However, the first region lacks the essential Glu active

site residue at position 241, and the second one lacks the essential Glu active site residue at

position 889. These domains are therefore predicted to be inactive.

細胞内局在 Cell membrane.

#### 画像

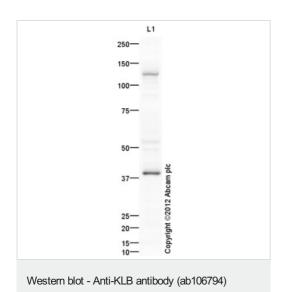

Anti-KLB antibody (ab106794) at 1  $\mu$ g/ml + HepG2 whole cell lysate at 10  $\mu$ g

#### **Secondary**

Goat Anti-Rabbit IgG H&L (HRP) preadsorbed (ab97080) at 1/5000 dilution

Developed using the ECL technique.


Performed under reducing conditions.

**Predicted band size:** 120 kDa **Observed band size:** 125 kDa

**Additional bands at:** 39 kDa, 54 kDa. We are unsure as to the identity of these extra bands.

Exposure time: 2 minutes

KLB contains a number of potential glycosylation sites (SwissProt) which may explain its migration at a higher molecular weight than predicted.

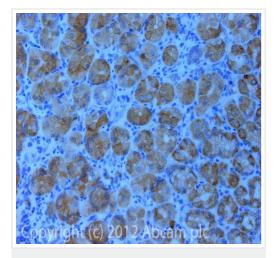

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-KLB antibody (ab106794)

IHC image of KLB staining in Human normal stomach formalin fixed paraffin embedded tissue section, performed on a Leica BondTM system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab106794, 5µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

# Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

### Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors