

Anti-Retinoic Acid Receptor beta antibody [EPR2017] ab124701

リコンビナント RabMAb


[2 Abreviews](#) [4 References](#) [画像数 7](#)

製品の概要

製品名	Anti-Retinoic Acid Receptor beta antibody [EPR2017]
製品の詳細	Rabbit monoclonal [EPR2017] to Retinoic Acid Receptor beta
由来種	Rabbit
アプリケーション	適用あり: WB, IHC-P, ICC/IF 適用なし: Flow Cyt or IP
種交差性	交差種: Human 非交差種: Mouse, Rat
免疫原	Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.
ポジティブ・コントロール	WB: Daudi, MCF7, MDA-MB-435, and U-87 MG lysates. IHC-P: human bladder carcinoma and kidney tissues. ICC/IF: HCT 116 cells.
特記事項	This product is a recombinant monoclonal antibody, which offers several advantages including: <ul style="list-style-type: none"> - High batch-to-batch consistency and reproducibility - Improved sensitivity and specificity - Long-term security of supply - Animal-free production For more information see here . Our RabMAb [®] technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to RabMAb[®] patents .

製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C. Stable for 12 months at -20°C.
解離定数 (K _D 値)	K _D = 1.13 x 10 ⁻¹¹ M

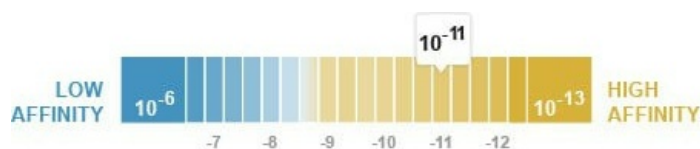

[Learn more about K_D](#)

バッファー	pH: 7.20 Preservative: 0.01% Sodium azide Constituents: 59% PBS, 40% Glycerol (glycerin, glycerine), 0.05% BSA
精製度	Protein A purified
ポリ/モノ	モノクローナル
クローン名	EPR2017
アイソタイプ	IgG

アプリケーション

The Abpromise guarantee **Abpromise保証は、** 次のテスト済みアプリケーションにおけるab124701の使用に適用されます
アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。


アプリケーション	Abreviews	特記事項
WB		1/1000 - 1/10000. Detects a band of approximately 55 kDa (predicted molecular weight: 50 kDa).
IHC-P		1/100 - 1/250. Perform heat mediated antigen retrieval before commencing with IHC staining protocol.
ICC/IF		1/100 - 1/250.

追加情報 Is unsuitable for Flow Cyt or IP.

ターゲット情報


機能	Receptor for retinoic acid. Retinoic acid receptors bind as heterodimers to their target response elements in response to their ligands, all-trans or 9-cis retinoic acid, and regulate gene expression in various biological processes. The RXR/RAR heterodimers bind to the retinoic acid response elements (RARE) composed of tandem 5'-AGGTCA-3' sites known as DR1-DR5. In the absence or presence of hormone ligand, acts mainly as an activator of gene expression due to weak binding to corepressors. In concert with RARG, required for skeletal growth, matrix homeostasis and growth plate function.
関連疾患	Microphthalmia, syndromic, 12
配列類似性	Belongs to the nuclear hormone receptor family. NR1 subfamily. Contains 1 nuclear receptor DNA-binding domain.
ドメイン	Composed of three domains: a modulating N-terminal domain, a DNA-binding domain and a C-terminal ligand-binding domain.
細胞内局在	Cytoplasm and Nucleus.

画像

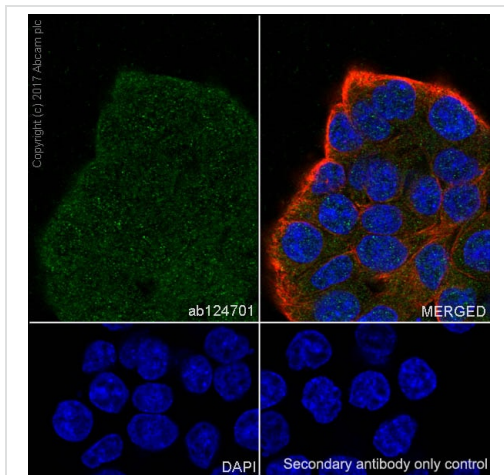


Immunocytochemistry/ Immunofluorescence - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

Immunocytochemistry/ Immunofluorescence analysis of HCT 116 (Human colorectal carcinoma epithelial cell) cells labeling Retinoic Acid Receptor beta with Purified ab124701 at 1:100 dilution (10 µg/ml). Cells were fixed in 4% Paraformaldehyde and permeabilized with 0.1% tritonX-100. Cells were counterstained with Ab195889 Anti-alpha Tubulin antibody [DM1A] - Microtubule Marker (Alexa Fluor® 594) 1:200 (2.5 µg/ml). Goat anti rabbit IgG (Alexa Fluor® 488, **ab150077**) was used as the secondary antibody at 1:1000 (2 µg/ml) dilution. DAPI (blue) was used as nuclear counterstain. PBS instead of the primary antibody was used as the secondary antibody only control.

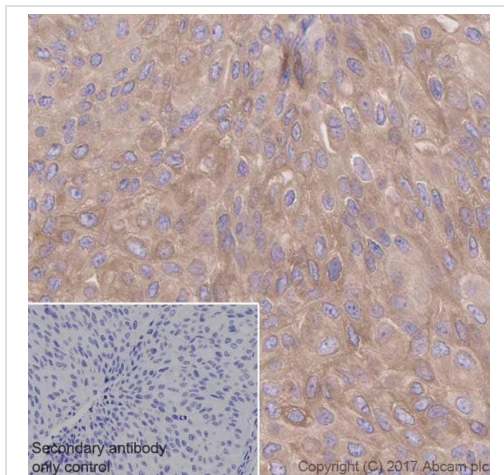

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) analysis of human bladder carcinoma tissue sections labeling Retinoic Acid Receptor beta with Purified ab124701 at 1:100 dilution (10.3 µg/ml). Perform heat mediated antigen retrieval using **ab93684** (Tris/EDTA buffer, pH 9.0). Purified)ImmunoHistoProbe one step HRP Polymer (ready to use) was used for detection. Negative control: PBS instead of the primary antibody. Hematoxylin was used as a counterstain.

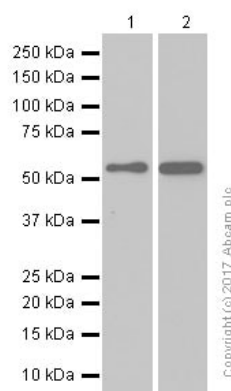

Western blot - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

All lanes : Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701) at 1/2000 dilution (Purified)

Lane 1 : Daudi (Human Burkitt's lymphoma lymphoblast) whole cell lysates

Lane 2 : U-87 MG (Human glioblastoma-astrocytoma epithelial cell) whole cell lysates


Lysates/proteins at 15 µg per lane.

Secondary

All lanes : Goat Anti-Rabbit IgG H&L (HRP) ([ab97051](#)) at 1/20000 dilution

Predicted band size: 50 kDa

Observed band size: 55 kDa

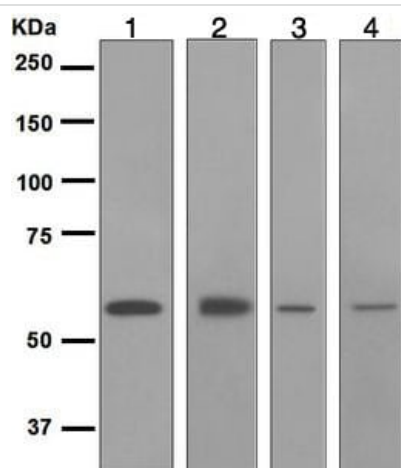

Western blot - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

All lanes : Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701) at 1/1000 dilution

Lane 1 : MCF7 cell lysates

Lane 2 : MDA-MB-435 cell lysates

Lane 3 : Daudi cell lysates

Lane 4 : U-87 MG cell lysates


Lysates/proteins at 10 µg per lane.

Secondary

Lane 1 : Goat anti-Rabbit HRP at 1/20000 dilution

Lanes 2-4 : Goat anti-Rabbit HRP at 1/2000 dilution


Predicted band size: 50 kDa

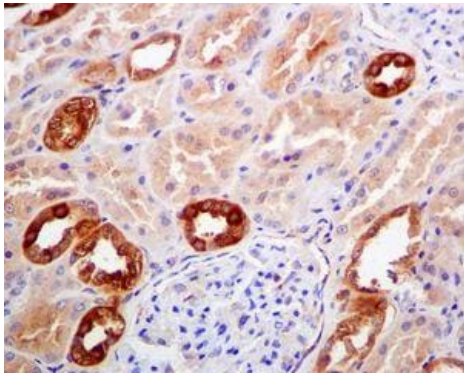

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

ab124701 at 1/100 dilution staining Retinoic Acid Receptor beta in paraffin-embedded Human kidney tissue by Immunohistochemistry.

Perform heat mediated antigen retrieval before commencing with IHC staining protocol.

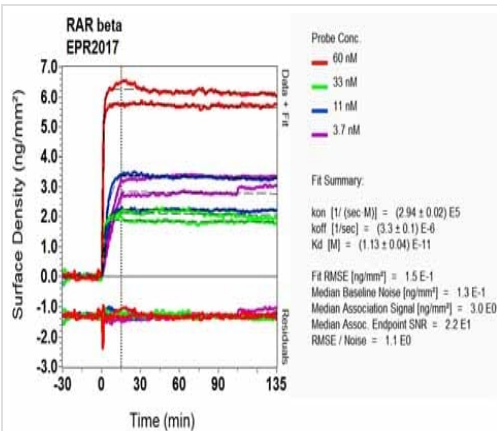


OIR-D Scanning - Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

Equilibrium dissociation constant (K_D)

Learn more about K_D

[Click here to learn more about \$K_D\$](#)

Why choose a recombinant antibody?

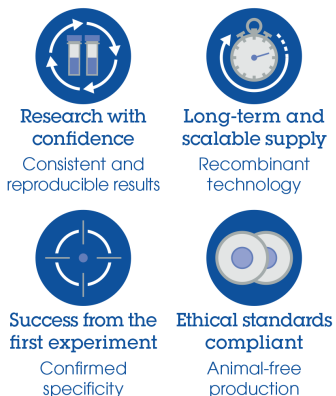

Anti-Retinoic Acid Receptor beta antibody [EPR2017] (ab124701)

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours

- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors