abcam

Product datasheet

Anti-VDAC1/Porin antibody ab34726

★★★★★ 6 Abreviews 59 References 画像数 1

製品の概要

製品名 Anti-VDAC1/Porin antibody

製品の詳細 Rabbit polyclonal to VDAC1/Porin

由来種 Rabbit

特異性 Cross reactivity will occur with all forms of VDACs including VDAC1, VDAC2 (4 isoforms) and

VDAC3(2 isoforms). This broad reactivity makes ab34726 useful as a loading control

(mitochondrial).

アプリケーション **適用あり**: IHC-P

種交差性 交差種: Human

交差が予測される動物種: Mouse, Rabbit, Cow, Pig, Zebrafish 🔷

免疫原 Synthetic peptide corresponding to Human VDAC1/Porin aa 185-197.

Sequence:

NDGTEFGGSIYQK

Run BLAST with
Run BLAST with

特記事項

The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Upon delivery aliquot and store at -20°C or -80°C. Avoid repeated freeze / thaw

cycles.

パッファー Preservative: 0.01% Sodium azide

Constituents: 0.42% Potassium phosphate, 0.87% Sodium chloride

精製度 Immunogen affinity purified

特記事項(精製) This antibody was affinity purified from monospecific antiserum.

1

ポリ/モノ ポリクローナル

アイソタイプ lgG

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab34726の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		Use a concentration of 5 µg/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

ターゲット情報

機能 Forms a channel through the mitochondrial outer membrane and also the plasma membrane. The

channel at the outer mitochondrial membrane allows diffusion of small hydrophilic molecules; in the plasma membrane it is involved in cell volume regulation and apoptosis. It adopts an open conformation at low or zero membrane potential and a closed conformation at potentials above 30-40 mV. The open state has a weak anion selectivity whereas the closed state is cation-selective. May participate in the formation of the permeability transition pore complex (PTPC)

responsible for the release of mitochondrial products that triggers apoptosis.

組織特異性 Heart, liver and skeletal muscle.


配列類似性 Belongs to the eukaryotic mitochondrial porin family.

ドメイン Consists mainly of a membrane-spanning beta-barrel formed by 19 beta-strands. The helical N-

terminus folds back into the pore opening and plays a role in voltage-gated channel activity.

細胞内局在 Mitochondrion outer membrane. Cell membrane.

画像

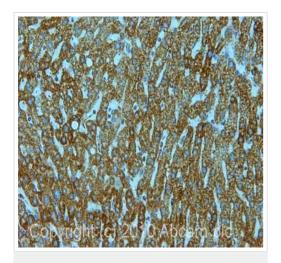

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-VDAC1/Porin antibody (ab34726)

IHC image of ab34726 staining in human liver formalin fixed paraffin embedded tissue section, performed on a Leica BondTM system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab34726, 5µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors