abcam

Product datasheet

Anti-Vitamin D Binding protein antibody ab65636

3 Abreviews 5 References 画像数 4

製品の概要

製品名 Anti-Vitamin D Binding protein antibody

製品の詳細 Rabbit polyclonal to Vitamin D Binding protein

由来種 Rabbit

アプリケーション 適用あり: IHC-P, ICC/IF, WB

種交差性 交差種: Human

交差が予測される動物種: Pig 🔷

免疫原 Synthetic peptide conjugated to KLH derived from within residues 450 to the C-terminus of

Human Vitamin D Binding Protein.lmmunogen の所有権に関して(Peptide available as ab71146.)

ポジティブ・コントロール WB: human liver and human ovary IHC-P: Human kidney FFPE tissue sections. ICC/IF: U937 cells

特記事項The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

バッファー pH: 7.40

Preservative: 0.02% Sodium azide

Constituent: PBS

Batches of this product that have a concentration < 1mg/ml may have BSA added as a stabilising

agent. If you would like information about the formulation of a specific lot, please contact our

scientific support team who will be happy to help.

精製度 Immunogen affinity purified

ポリ/モノ ポリクローナル

1

アプリケーション

The Abpromise guarantee Abpromise保証は、次のテスト済みアプリケーションにおけるab65636の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

| アプリケーション | Abreviews | 特記事項 |
|----------|-----------|--|
| IHC-P | | Use a concentration of 1 µg/ml. |
| ICC/IF | | Use a concentration of 5 µg/ml. |
| WB | | Use a concentration of 1 µg/ml. Detects a band of approximately 53 kDa (predicted molecular weight: 53 kDa). |

ターゲット情報

機能 Multifunctional protein found in plasma, ascitic fluid, cerebrospinal fluid, and urine and on the

surface of many cell types. In plasma, it carries the vitamin D sterols and prevents polymerization of actin by binding its monomers. DBP associates with membrane-bound immunoglobulin on the


surface of B-lymphocytes and with IgG Fc receptor on the membranes of T-lymphocytes.

配列類似性 Belongs to the ALB/AFP/VDB family.

Contains 3 albumin domains.

細胞内局在 Secreted.

画像

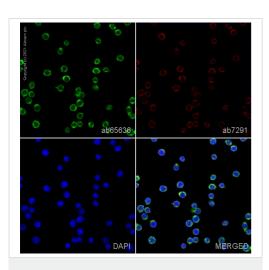


Immunocytochemistry/ Immunofluorescence - Anti-Vitamin D Binding protein antibody (ab65636) ab65636 staining Vitamin D binding protein in U937 cells. The cells were fixed with 100% methanol (5 min), permeabilized with 0.1% PBS-Tween for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at 4°C with ab65636 at 5µg/ml and ab7291, Mouse monoclonal [DM1A] to alpha Tubulin - Loading Control. Cells were then incubated with ab150081, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 488), preadsorbed at 1/1000 dilution (shown in green) and ab150080, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 594) at 1/1000 dilution (shown in pseudocolour red). Nuclear DNA was labelled with DAPI (shown in blue).

Image was acquired with a confocal microscope (Leica-Microsystems TCS SP8) and a single confocal section is shown.



Western blot - Anti-Vitamin D Binding protein antibody (ab65636)

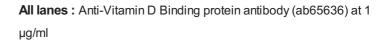

Lane 1 : Human liver tissue lysate - total protein (ab29889)

Lane 2: Human ovary tissue lysate - total protein (ab30222)

Lysates/proteins at 10 µg per lane.


Secondary

All lanes : Goat polyclonal to Rabbit lgG - H&L - Pre-Adsorbed (HRP) at 1/3000 dilution

Performed under reducing conditions.


Predicted band size: 53 kDa Observed band size: 53 kDa

Exposure time: 1 minute

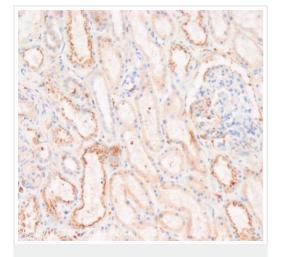

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Vitamin D Binding protein antibody (ab65636)

IHC image of Vitamin D Binding Protein staining in human kidney FFPE section, performed on a BondTM system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab65636, 1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX

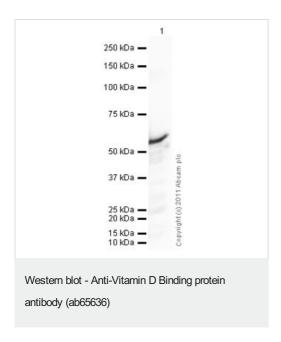

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors