abcam

Product datasheet

Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647) ab150075

★★★★★ 5 Abreviews 450 References 画像数 5

製品の概要

製品名 Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647)

由来種Donkeyターゲット生物種Rabbit

アプリケーション 適用あり: ICC/IF, ELISA, IHC-P, Flow Cyt, IHC-Fr

免疫原 The details of the immunogen for this antibody are not available.

標識 Alexa Fluor® 647. Ex: 652nm, Em: 668nm

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Stable for 12 months at -20°C. Store In the Dark.

ארעדעד Preservative: 0.02% Sodium azide

Constituents: 23% Glycerol (glycerin, glycerine), PBS, 1% BSA

精製度 Immunogen affinity purified

特記事項(精製) This antibody was isolated by affinity chromatography using antigen coupled to agarose beads.

ポリ/モノ ポリクローナル

アイソタイプ lgG

特記事項 Alexa Fluor[®] is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific

Company. The Alexa Fluor[®] dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor[®] dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor[®] dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact

Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or

1


outlicensing@thermofisher.com.

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab150075の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
ICC/IF	**** <u>(1)</u>	1/200 - 1/1000.
ELISA		Use at an assay dependent concentration.
IHC-P		Use at an assay dependent concentration.
Flow Cyt		1/2000 - 1/4000. ab199093 - Rabbit monoclonal lgG (Alexa Fluor® 647), is suitable for use as an isotype control to complement this secondary antibody.
IHC-Fr	★★★★★ (2)	Use at an assay dependent concentration.


画像

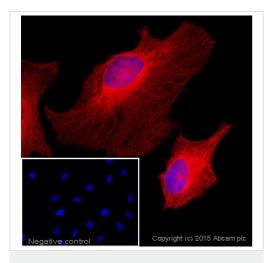


Immunocytochemistry/ Immunofluorescence Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647)
(ab150075)

ICC/IF image of <u>ab6046</u> stained HeLa cells. The cells were 100% formaldehyde fixed (10 min), permeabilized with 0.1% Triton X-100 for 5 minutes and then incubated in 1%BSA / 10% normal donkey serum / 0.3M glycine in 0.1% PBS-Tween for 1h to block non-specific protein-protein interactions. The cells were then incubated with the primary antibody (<u>ab6046</u>, 5 μ g/ml) overnight at +4°C. The secondary antibody (red) was ab150075 Alexa Fluor® 647 donkey anti-rabbit lgG (H+L) used at 1 μ g/ml for 1h.DAPI was used to stain the cell nuclei (blue) at a concentration of 1.43 μ M.

The negative control (inset) is a secondary-only assay to demonstrate low non-specific binding of the secondary antibody.

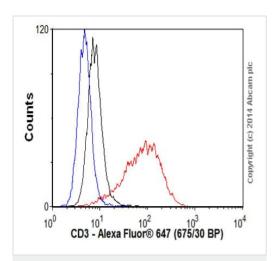

Flow Cytometry - Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647) (ab150075)

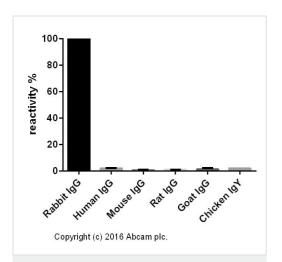


ELISA - Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647) (ab150075)

Overlay histogram showing Jurkat cells stained with <u>ab16669</u> (red line). The cells were fixed with 4% paraformaldehyde (10 min) and then permeabilized with 0.1% PBS-Tween for 20 min. The cells were then incubated in 1x PBS / 10% normal donkey serum / 0.3M glycine to block non-specific protein-protein interactions followed by the antibody (<u>ab16669</u>, 1/1000 dilution) for 30 min at 22°C. The secondary antibody Donkey anti-rabbit lgG H&L (Alexa Fluor® 647) (ab150075) was used at 1/2000 dilution for 30 min at 22°C. Isotype control antibody (black line) was rabbit lgG (monclonal) (0.1µg/1x10⁶ cells) used under the same conditions. Unlabelled sample (blue line) was also used as a control. Acquisition of >5,000 events were collected using a solid-state 25mW red diode laser (635nm) and 675/30 bandpass filter.

Cross-reactivity of the polyclonal secondary antibody <u>ab182020</u> was tested using a sandwich ELISA approach. The wells were coated with the indicated IgG standards at 1 μ g/ml (50 μ l/well) and incubated overnight at 4°C, followed by a 5% BSA blocking step for 2h at RT. <u>ab182020</u> was then added starting at 1 μ g/ml and gradually diluted 1/4 (50 μ l/well), followed by incubation for 2h. For the detection Goat anti-Donkey IgG H&L (HRP) (<u>ab6988</u>) was used at 1/20,000 dilution (50 μ l/well), followed by incubation for 1h at RT.

For the batch tested, <u>ab182020</u> showed a cross-reactivity below 2% towards human IgG, mouse IgG, rat IgG, goat IgG and chicken IgY.


This data was developed using the unconjugated antibody (ab182020).

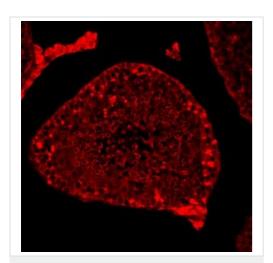

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647) (ab150075)

This image is courtesy of an anonymous Abreview.

<u>ab62050</u> staining prostaglandin E synthase in mouse testis tissue sections by Immunohistochemistry (IHC-P - paraformaldehydefixed, paraffin-embedded sections). Tissue was fixed with paraformaldehyde, permeabilized with Triton X-100, and blocked with 10% serum for 30 minutes at 24°C; antigen retrieval was by heat mediation in a citrate buffer. Samples were incubated with primary antibody (1/50) for 20 hours at 4°C. An Alexa Fluor® 647-conjugated donkey anti-rabbit IgG H&L (ab150075) (1/200) was used as the secondary antibody.

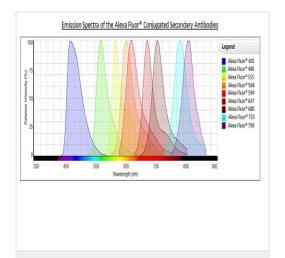

Alexa Fluor® - Donkey Anti-Rabbit IgG H&L (Alexa Fluor® 647) (ab150075)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- · Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors