# abcam

# Product datasheet

# GsMTx-4, mechanosensitive and stretch-activated ion channel inhibitor ab141871

★★★★★ 1 Abreviews 29 References 画像数 2

#### 製品の概要

製品名 GsMTx-4, mechanosensitive and stretch-activated ion channel inhibitor

製品の詳細 Mechanosensitive and stretch-activated ion channel inhibitor

**CAS 番号** 1209500-46-8

構造式 Gly-Cys-Leu-Glu-Phe-Trp-Trp-Lys-Cys-Asn-Pro-Asn-Asp-

Asp-Lys-Cys-Cys-Arg-Pro-Lys-Leu-Lys-Cys-Ser-Lys-Leu-

Phe-Lys-Leu-Cys-Asn-Phe-Ser-Phe-NH<sub>2</sub>

(Disulfide bridges: 2-17, 9-23 and 16-30)

#### 製品の特性

分子量 4095.85

分子式 C<sub>185</sub>H<sub>273</sub>N<sub>49</sub>O<sub>45</sub>S<sub>6</sub>

配列 GCLEFWWKCNPNDDKCCRPKLKCSKLFKLCNFSF (Modifications: Phe-34 = C-terminal

amide, Disulfide bridges: 2-17, 9-23 and 16-30)

**PubChem 登録番号** 90488987

保存方法 Store at -20°C. Store under desiccating conditions. The product can be stored for up to 12

months.

溶解性 Soluble in water to 1mg/ml

使用に関する注意 Wherever possible, you should prepare and use solutions on the same day. However, if you need

to make up stock solutions in advance, we recommend that you store the solution as aliquots in tightly sealed vials at -20°C. Generally, these will be useable for up to one month. Before use, and

prior to opening the vial we recommend that you allow your product to equilibrate to room

temperature for at least 1 hour.

Need more advice on solubility, usage and handling? Please visit our frequently asked

questions (FAQ) page for more details.


由来 Grammostola rosea


画像

Gly-Cys-Leu-Glu-Phe-Trp-Trp-Lys-Cys-Asn-Pro-Asn-Asp-Asp-Lys-Cys-Cys-Arg-Pro-Lys-Leu-Lys-Cys-Ser-Lys-Leu-Phe-Lys-Leu-Cys-Asn-Phe-Ser-Phe-NH<sub>2</sub>

(Disulfide bridges: 2-17, 9-23 and 16-30)

Chemical Structure - GsMTx-4, mechanosensitive and stretch-activated ion channel inhibitor (ab141871)

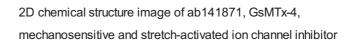

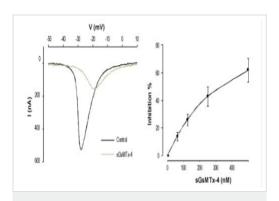

Functional Studies - GsMTx-4, mechanosensitive and stretch-activated ion channel inhibitor (ab141871) GsMTx-4 (<u>ab141870</u>) inhibits Na $_{
m V}$ 1.7 channel currents expressed in *Xenopus* oocytes. Na $_{
m V}$ 1.7 currents were elicited by 100 ms voltage ramp from a holding potential of -100 mV to +30 mV, applied every 10 seconds using whole-cell voltage clamp technique. Left: Superimposed traces of Na $_{
m V}$ 1.7 currents before (black) and during (green) application of 500 nM GsMTx-4 (<u>ab141870</u>). Right: GsMTx-4 dose response inhibition of Na $_{
m V}$ 1.7 currents.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES, NOT FOR USE IN HUMANS"

### Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- · We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

## Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors
- Abcam biochemicals are novel compounds and we have not tested their biological activity in house. Please use the literature to identify how to use these products effectively. If you require further assistance please contact the scientific support team