abcam

Product datasheet

Proteinase K Antigen Retrieval Solution ab64220

★★★★★ 1 Abreviews 24 References 画像数 1

製品の概要

製品名

アプリケーション

特記事項

Proteinase K Antigen Retrieval Solution

適用あり: IHC-P, In situ hybridization

ab64220 is a single-component, ready-to-use proteinase K solution for use in protease-induced antigen retrieval during IHC (PIER) and pre-treatment for in situ hybridization.

Staining Protocol:

IHC procedure: Incubate tissue section for 5 to 7 minutes at room temperature; Rinse and proceed with immunostaining.

General ISH procedure: Incubate tissue section for 15 minutes at 37°C; Rinse and proceed with hybridization.

Other products for IHC

Other enzymes for antigen retrieval include: <u>Pepsin Solution ab64201</u> and <u>Trypsin Antigen</u> Retrieval Kit ab259256.

Find more kits and reagents for antigen retrieval, blocking, signal amplification, visualization, counterstaining, and mounting in the **IHC kits and reagents guide**.

製品の特性

製品の状態 Liquid

保存方法 Store at -20°C. Stable for 12 months at -20°C

バッファー pH: 7

Preservative: 0.1% Sodium azide

Constituents: Tris buffer, 10% Glycerol (glycerin, glycerine)

関連性 Proteinase K is an endolytic protease that cleaves peptide bonds preferentially next to the

carboxyl group of N-substituted hydrophobic aliphatic and aromatic amino acids. Proteinase K is classified as a serine protease and exhibits broad substrate specificity. It can hydrolyze keratin; hence the name. Proteinase K is used for the digestion of proteins on tissue sections as a pretreatment to aid in the preparation of mRNA and its detection by in situ hybridisation methods


using oligonucleotide probes.

1

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab64220の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		Use at an assay dependent dilution.
In situ hybridization		Use at an assay dependent dilution.

画像

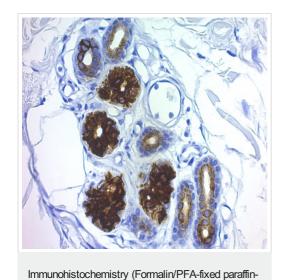

ab64220 Proteinase K used with aa anti-Pan CK antibody on Human skin

embedded sections) - Proteinase K (ab64220)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

Guarantee only valid for products bought direct from Abcam or one of our authorized distributors				
		3		