abcam

Product datasheet

Anti-3-Nitrotyrosine antibody [7A12AF6] ab110282

1 Abreviews 26 References 画像数 5

製品の概要

製品名 Anti-3-Nitrotyrosine antibody [7A12AF6]

製品の詳細 Mouse monoclonal [7A12AF6] to 3-Nitrotyrosine

由来種 Mouse

特異性 ab110282 was developed to recognize only protein-bound nitrotyrosine and so is a sensitive tool

for measuring protein-specific modifications from oxidative stress.

アプリケーション 適用あり: ICC/IF, Flow Cyt, WB, IP, In-Cell ELISA

種交差性 交差種: Species independent

免疫原 Chemical/ Small Molecule. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール Flow Cyt: HL-60 cells treated with 2 mM peroxynitrite ICC/IF: HeLa cells and Human fibroblast

cells treated with 1 mM peroxynitrite WB: nitrated Bovine heart mitochondria; nitrated tyrosine

特記事項

This antibody clone is manufactured by Abcam. If you require a custom buffer formulation or

conjugation for your experiments, please contact **orders@abcam.com**.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

Product was previously marketed under the MitoSciences sub-brand.

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C. Do Not Freeze.

バッファー pH: 7.5

Preservative: 0.02% Sodium azide Constituent: HEPES buffered saline

精製度 Proprietary Purification

特記事項(精製) The antibody was produced in vitro using hybridomas grown in serum-free medium, and then

1

purified by biochemical fractionation. Purity >95% by SDS-PAGE.

ポリ/モノ モノクローナル

クローン名 7A12AF6

アイソタイプ lgG2b

軽鎖の種類 kappa

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab110282の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。


アプリケーション	Abreviews	特記事項
ICC/IF		Use a concentration of 1 µg/ml.
Flow Cyt		Use a concentration of 1 μ g/ml. <u>ab170192</u> - Mouse monoclonal lgG2b, is suitable for use as an isotype control with this antibody.
WB		Use a concentration of 1 µg/ml.
IP		Use at an assay dependent concentration.
In-Cell ELISA		Use a concentration of 4 µg/ml. (0.4 µg/well).

ターゲット情報

関連性

Protein tyrosine nitration results in a post-translational modification that is increasingly receiving attention as an important component of nitric oxide signaling. While multiple nonenzymatic mechanisms are known to be capable of producing nitrated tyrosine residues, most tyrosine nitration events involve catalysis by metalloproteins such as myeloperoxidase, eosinophilperoxidase, myoglobin, the cytochrome P-450s, superoxide dismutase and prostacyclin synthase. Various studies have shown that protein tyrosinenitration is limited to specific proteins and that the process is selective. For example, exposure of human surfactant protein A, SP-A, to oxygen-nitrogen intermediates generated by activated alveolar macrophages resulted in specific nitration of SP-A at tyrosines 164 and 166, while addition of 1.2 mMCO 2 resulted in additional nitration at tyrosine 161. The presence of nitrotyrosine-containing proteins has shown high correlation to disease states such as atherosclerosis, Alzheimer's disease, Parkinson's disease and amyotrophic lateral sclerosis.

画像

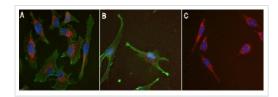


Immunocytochemistry/ Immunofluorescence - Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282)

Immunocytochemistry image of ab110282 stained Human HeLa cells (A) and fibroblast cells (B, C).


Cells grown on slides were paraformaldehyde fixed (4%, 20 min) and Triton X-100 permeabilized (0.1%, 15 min). Slides were treated with/without 1 mM peroxynitrite to modify exposed tyrosines to 3-nitrotyrosine. Slides were blocked and incubated with tab110282 at 2 μ g/ml overnight at 4°C. The secondary antibody (green) was Alexa Fluor® 488 Goat anti-Mouse lgG (H+L) used at a 1/1000 dilution for 1 hour. 10% Goat serum was used as the blocking agent for all blocking steps. For reference, the mitochondria (red) were identified by HSP60 / Alexa Fluor® 594 and DAPI was used to stain the cell nuclei (blue).

HeLa cells (A) and fibroblast cells (B) show surface modification of tyrosine to 3-nitrotyrosine after exposure to peroxynitrite. While (C) unexposed fibroblast cells show no modification.

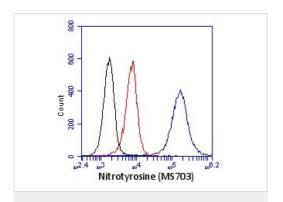

Flow Cytometry - Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282)

HL-60 cells were stained with 1 μ g/ml ab110282 following treatment with 2 mM peroxynitrite (blue) or vehicle control (red). No primary antibody control is shown in black. Peroxynitrite modifies tyrosine residues to 3-nitrotyrosine.

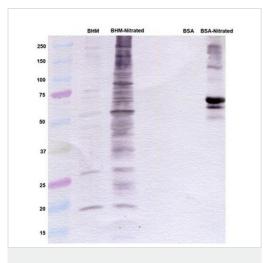

Western blot - Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282)


All lanes : Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282) at 1 µg/ml

Lane 1: Bovine heart mitochondria


Lane 2: Bovine heart mitochondria - nitrated

Lane 3: BSA

Lane 4: BSA - nitrated

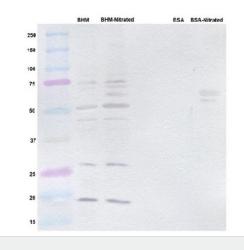

Western blot - Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282)

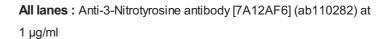


Lane 1: Bovine heart mitochondria

Lane 2: Bovine heart mitochondria - nitrated

Lane 3: BSA

Lane 4: BSA - nitrated

Prior to running the samples, the membrane was treated with sodium dithionite to reduce nitrotyrosine to aminotyrosine.

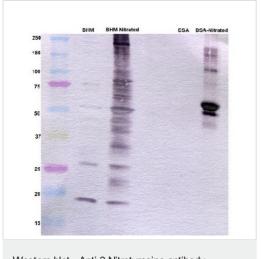

Western blot - Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282)

All lanes : Anti-3-Nitrotyrosine antibody [7A12AF6] (ab110282) at 1 μ g/ml

Lane 1: Bovine heart mitochondria

Lane 2: Bovine heart mitochondria - nitrated

Lane 3: BSA

Lane 4: BSA - nitrated

In this experiment, the antibody was first blocked with free nitrotyrosine before being used to blot the membrane.

The figure shows that ab110282's binding capacity was not inhibited by the free nitrotyrosine, and so only binds to the protein-bound form.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors