# abcam

### Product datasheet

## Alexa Fluor® 488 Anti-BMP4 antibody [EPR6211] ab200794

יעלאעבע RabMAb

#### 画像数3

#### 製品の概要

特記事項

製品名 Alexa Fluor® 488 Anti-BMP4 antibody [EPR6211]

製品の詳細 Alexa Fluor® 488 Rabbit monoclonal [EPR6211] to BMP4

由来種 Rabbit

標識 Alexa Fluor® 488. Ex: 495nm, Em: 519nm

アプリケーション 適用あり: Flow Cyt (Intra), ICC/IF

種交差性 交差種: Human

免疫原 Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール ICC/IF: HeLa cells. Flow Cyt (intra): HeLa cells.

This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility

- Improved sensitivity and specificity

- Long-term security of supply

- Animal-free production

For more information see here.

Our RabMAb® technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb patents**.


Alexa Fluor<sup>®</sup> is a registered trademark of Molecular Probes, Inc., a Thermo Fisher Scientific Company. The Alexa Fluor® dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor® dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor<sup>®</sup> dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@thermofisher.com.

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

解離定数( $K_D$ 値)  $K_D = 1.94 \times 10^{-10} M$ 

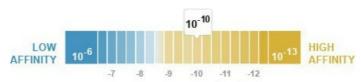

Learn more about K<sub>D</sub>

**バッファー** pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 30% Glycerol (glycerin, glycerine), 1% BSA, PBS

精製度 Protein A purified

**ポリ/モノ** モノクローナル

**クローン名** EPR6211

アイソタイプ IgG

#### アプリケーション

**The Abpromise guarantee** <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab200794の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
Flow Cyt (Intra)		1/50.
ICC/IF		1/100.

#### ターゲット情報

機能 Induces cartilage and bone formation. Also act in mesoderm induction, tooth development, limb

formation and fracture repair. Acts in concert with PTHLH/PTHRP to stimulate ductal outgrowth

during embryonic mammary development and to inhibit hair follicle induction.

組織特異性 Expressed in the lung and lower levels seen in the kidney. Present also in normal and neoplastic

prostate tissues, and prostate cancer cell lines.

**関連疾患** Defects in BMP4 are the cause of microphthalmia syndromic type 6 (MCOPS6) [MIM:607932];

also known as microphthalmia and pituitary anomalies or microphthalmia with brain and digit developmental anomalies. Microphthalmia is a clinically heterogeneous disorder of eye formation,

ranging from small size of a single eye to complete bilateral absence of ocular tissues (anophthalmia). In many cases, microphthalmia/anophthalmia occurs in association with

syndromes that include non-ocular abnormalities. MCOPS6 is characterized by

microphthalmia/anophthalmia associated with facial, genital, skeletal, neurologic and endocrine

anomalies.

Defects in BMP4 are the cause of non-syndromic orofacial cleft type 11 (OFC11) [MIM:600625].

Non-syndromic orofacial cleft is a common birth defect consisting of cleft lips with or without cleft palate. Cleft lips are associated with cleft palate in two-third of cases. A cleft lip can occur on one or both sides and range in severity from a simple notch in the upper lip to a complete opening in the lip extending into the floor of the nostril and involving the upper gum. OFC11 is an unusual anomaly consisting of a paramedian scar of the upper lip with an appearance suggesting that a typical cleft lip was corrected in utero.


配列類似性

細胞内局在

Belongs to the TGF-beta family.

Secreted > extracellular space > extracellular matrix.


#### 画像

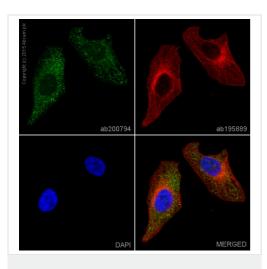

Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 488 Anti-BMP4 antibody [EPR6211] (ab200794)


ab200794 staining BMP4 in HeLa cells. The cells were fixed with 4% formaldehyde (10min), permeabilized with 0.1% Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at +4°C with ab200794 at 1/100 dilution (shown in green) and <u>ab195889</u>, Mouse monoclonal to alpha Tubulin (Alexa Fluor<sup>®</sup> 594), at  $2\mu g/ml$  (shown in red). Nuclear DNA was labelled with DAPI (shown in blue).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

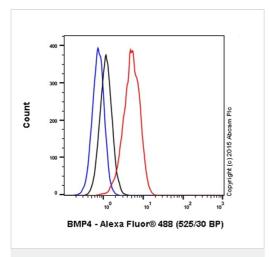

Flow Cytometry (Intracellular) - Alexa Fluor® 488 Anti-BMP4 antibody [EPR6211] (ab200794)

Overlay histogram showing HeLa cells stained with ab200794 (red line). The cells were fixed with 4% formaldehyde (10 min) and then permeabilized with 0.1% PBS-Tween for 20 min. The cells were then incubated in 1x PBS / 10% normal goat serum / 0.3M glycine to block non-specific protein-protein interactions followed by the antibody (ab200794, 1/50 dilution) for 30 min at 22°C. Isotype control antibody (black line) was rabbit monoclonal IgG [EPR25A] Alexa Fluor® 488 (ab199091) used at the same concentration and conditions as the primary antibody. Unlabelled sample (blue line) was also used as a control. Acquisition of >5,000 events were collected using a 20mW Argon ion laser (488nm) and 525/30 bandpass filter. This antibody gave a positive signal in HeLa cells fixed with 80% methanol (5 min)/permeabilized with 0.1% PBS-Tween for 20 min used under the same conditions.

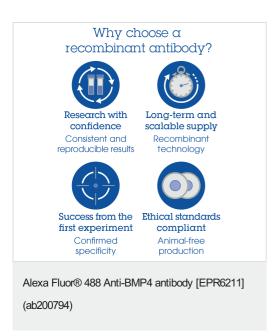

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

#### Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

#### Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors