# abcam

## Product datasheet

## Alexa Fluor® 488 Anti-CD4 antibody [EPR19514] ab277275

リコンピナント

★★★★★ 1 Abreviews 画像数 2

#### 製品の概要

特記事項

製品名 Alexa Fluor® 488 Anti-CD4 antibody [EPR19514]

製品の詳細 Alexa Fluor® 488 Rabbit monoclonal [EPR19514] to CD4

由来種 Rabbit

標識 Alexa Fluor® 488. Ex: 495nm, Em: 519nm

アプリケーション **適用あり**: IHC-P

種交差性 交差種: Mouse

免疫原 Recombinant fragment. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール IHC-P: Mouse normal spleen tissue.

This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility
- Improved sensitivity and specificity
- Long-term security of supply
- Animal-free production

For more information see here.

Our RabMAb<sup>®</sup> technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb**<sup>®</sup> **patents**.

Alexa Fluor<sup>®</sup> is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor<sup>®</sup> dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor<sup>®</sup> dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor<sup>®</sup> dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or **outlicensing@thermofisher.com**.

#### 製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long

term. Avoid freeze / thaw cycle. Store In the Dark.

**バッファー** pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 68% PBS, 30% Glycerol (glycerin, glycerine), 1% BSA

精製度 Protein A purified

**ポリ/モノ** モノクローナル **クローン名** EPR19514

アイソタイプ IgG

#### アプリケーション

**The Abpromise guarantee** <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab277275の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P	<b>★</b> ★ ★ ★ ★ <b>(1)</b>	1/100. Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

## ターゲット情報

機能 Accessory protein for MHC class-II antigen/T-cell receptor interaction. May regulate T-cell

activation. Induces the aggregation of lipid rafts.

配列類似性 Contains 3 lg-like C2-type (immunoglobulin-like) domains.

Contains 1 lg-like V-type (immunoglobulin-like) domain.

翻訳後修飾 Palmitoylation and association with LCK contribute to the enrichment of CD4 in lipid rafts.


**細胞内局在** Cell membrane. Localizes to lipid rafts. Removed from plasma membrane by HIV-1 Nef protein

that increases clathrin-dependent endocytosis of this antigen to target it to lysosomal degradation.

Cell surface expression is also down-modulated by HIV-1 Envelope polyprotein gp160 that

interacts with, and sequesters CD4 in the endoplasmic reticulum.

#### 画像

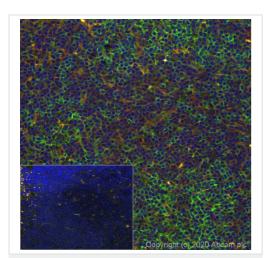


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Alexa Fluor® 488 Anti-CD4 antibody [EPR19514] (ab277275)

Immunofluorescence analysis of CD4 staining in a section of formalin-fixed paraffin-embedded mouse normal spleen.

Performed on a Leica BOND. The section was pre-treated using heat mediated antigen retrieval with EDTA (pH 9.0) using retrieval settings of 110°C for 20 minutes. The section was then incubated at room temperature for 1 hour with ab277275 at 1/100 dilution (shown in green). Nuclear DNA was labeled with DAPI (shown in blue). The section was then mounted using Dako Fluorescence Mounting Medium<sup>®</sup>.

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

For other IHC staining systems (automated and non-automated), customers should optimize variable parameters such as antigen retrieval conditions, antibody concentrations and incubation times.

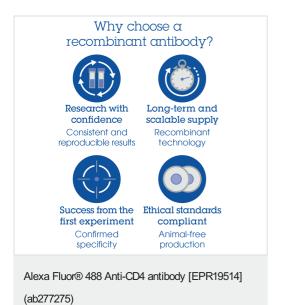

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

## Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery

- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

## Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors