abcam

Product datasheet

Alexa Fluor® 647 Anti-BrdU antibody [BU1/75 (ICR1)] ab220075

4 References 画像数 1

製品の概要

特記事項

製品名 Alexa Fluor® 647 Anti-BrdU antibody [BU1/75 (ICR1)]

製品の詳細 Alexa Fluor® 647 Rat monoclonal [BU1/75 (ICR1)] to BrdU

由来種 Rat

標識 Alexa Fluor® 647. Ex: 652nm, Em: 668nm

特異性
This antibody reacts with BrdU in single stranded DNA, BrdU attached to a protein carrier or free BrdU. It detects nucleated cells in S-Phase which have had BrdU incorporated into their DNA.
Also reacts with chlorodeoxyuridine but with reduced staining. The antibody does not react with

thymidine. It has been reported in the literature that this antibody clone cross-reacts with Edu

(PMID: 23272138) and some customers reported that it cross reacts with IdU.

アプリケーション 適用あり: ICC/IF

種交差性 交差種: Species independent

免疫原 The details of the immunogen for this antibody are not available.

ポジティブ・コントロール ICC/IF: HeLa cells treated with BrdU (10uM for 24 hours)

The antibody recognises single stranded DNA so the DNA needs to be unraveled first. This can be done with DNAse, although this doesn't give the best results. Depending on the assay, acid denaturation with 2M HCL or heat denaturation are the most successful. Please note this step is critical in any assay with this antibody and is the area that should be modified to optimise results. Detailed BrdU protocol is available in "Neuroscience protocols" on our "Protocol and

troubleshooting tips" webpage.

Alexa Fluor[®] is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor[®] dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor[®] dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor[®] dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or

outlicensing@thermofisher.com.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

バッファー pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 30% Glycerol (glycerin, glycerine), 1% BSA, PBS

精製度 Immunogen affinity purified

ポリ/モノ モノクローナル **DU1/75 (ICR1)**

アイソタイプ lgG2a

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab220075の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。


アプリケーション	Abreviews	特記事項
ICC/IF		1/100. This product gave a positive signal in HeLa cells treated with BrdU (10uM for 24 hours) fixed with 4% formaldehyde (10 min) and 100% methanol (5 min)

ターゲット情報

関連性

The immunocytochemical detection of bromodeoxyuridine (BrdU) incorporated into DNA is a powerful tool to study the cytokinetics of normal and neoplastic cells. In vitro or in vivo labeling of tumor cells with the thymidine analogue BrdU and the subsequent detection of incorporated BrdU with specific anti-BrdU monoclonal antibodies is an accurate and comprehensive method to quantitate the degree of DNA-synthesis. BrdU is incorporated into the newly synthezised DNA of S-phase cells may provide an estimate for the fraction of cells in S-phase. Also dynamic proliferative information such as the S-phase transit rate and the potential doubling time can be obtained, by means of bivariate BrdU/DNA flow cytometric analysis.

細胞内局在 Nuclear

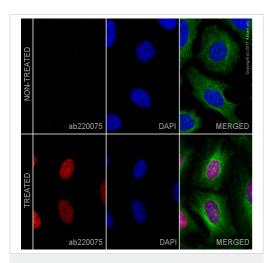

Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 647 Anti-BrdU antibody [BU1/75 (ICR1)] (ab220075)

ab220075 staining BrdU in HeLa cells. Untreated and BrdU treated (ab142567, 10uM for 24 hours) cells were fixed with 100% methanol (5 min) and then subjected to acid hydrolysis using 2M HCL in 0.1% PBS-Tween for 30 minutes at room temperature to denature the DNA. They were then permeabilized with 0.1% Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at +4°C with ab220075 at 1/100 dilution (shown in red) and ab195887, Mouse monoclonal to alpha Tubulin (Alexa Fluor® 488), at 1/250 dilution (shown in green). Nuclear DNA was labelled with DAPI (shown in blue).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

This product also gave a positive signal under the same testing conditions in HeLa cells fixed with 4% formaldehyde (10 min).

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

Guarantee only valid for products bought direct from Abcam or one of our authorized distributors