# abcam

### **Product datasheet**

# Alexa Fluor® 647 Anti-Caldesmon/CDM antibody [E89] ab208117

אילשעבע RabMAb

画像数 2

製品の概要		
製品名	Alexa Fluor® 647 Anti-Caldesmon/CDM antibody [E89]	
製品の詳細	Alexa Fluor® 647 Rabbit monoclonal [E89] to Caldesmon/CDM	
由来種	Rabbit	
標識	Alexa Fluor® 647. Ex: 652nm, Em: 668nm	
アプリケーション	<b>適用あり:</b> ICC/IF	
種交差性	交差種: Human	
免疫原	Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.	
ポジティブ・コントロール	ICC/IF: HeLa cells.	
特記事項	Our RabMAb <sup>®</sup> technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to <b><u>RabMAb<sup>®</sup> patents</u></b> .	
	Alexa Fluor <sup>®</sup> is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor <sup>®</sup> dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor <sup>®</sup> dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor <sup>®</sup> dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or <b>outlicensing@thermofisher.com</b> .	

製品の特性

製品の状態

Liquid

保存方法

Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

	Avoid freeze / thaw cycle. Store In the Dark.
バッファー	pH: 7.40 Preservative: 0.02% Sodium azide Constituents: 30% Glycerol (glycerin, glycerine), 1% BSA, PBS
精製度	Protein A purified
ポリ/モノ	モノクローナル
クローン名	E89
アイソタイプ	lgG

アプリケーション


## The Abpromise guarantee Abpromise保証は、次のテスト済みアプリケーションにおけるab208117の使用に適用されます

#### アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
ICC/IF		1/100. This product gave a positive signal in HeLa cells fixed with 4% formaldehyde (10 min).


#### ターゲット情報

機能	Actin- and myosin-binding protein implicated in the regulation of actomyosin interactions in smooth muscle and nonmuscle cells (could act as a bridge between myosin and actin filaments). Stimulates actin binding of tropomyosin which increases the stabilization of actin filament structure. In muscle tissues, inhibits the actomyosin ATPase by binding to F-actin. This inhibition is attenuated by calcium-calmodulin and is potentiated by tropomyosin. Interacts with actin, myosin, two molecules of tropomyosin and with calmodulin. Also play an essential role during cellular mitosis and receptor capping.
組織特異性	High-molecular-weight caldesmon (isoform 1) is predominantly expressed in smooth muscles, whereas low-molecular-weight caldesmon (isoforms 2, 3, 4 and 5) are widely distributed in non-muscle tissues and cells. Not expressed in skeletal muscle or heart.
配列類似性	Belongs to the caldesmon family.
ドメイン	The N-terminal part seems to be a myosin/calmodulin-binding domain, and the C-terminal a tropomyosin/actin/calmodulin-binding domain. These two domains are separated by a central helical region in the smooth-muscle form.
翻訳後修飾	In non-muscle cells, phosphorylation by CDK1 during mitosis causes caldesmon to dissociate from microfilaments. Phosphorylation reduces caldesmon binding to actin, myosin, and calmodulin as well as its inhibition of actomyosin ATPase activity. Phosphorylation also occurs in both quiescent and dividing smooth muscle cells with similar effects on the interaction with actin and calmodulin and on microfilaments reorganization.
細胞内局在	Cytoplasm > cytoskeleton. Cytoplasm > myofibril. On thin filaments in smooth muscle and on stress fibers in fibroblasts (nonmuscle).

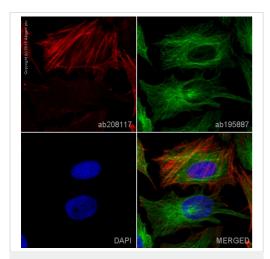

Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 647 Anti-Caldesmon/CDM antibody [E89] (ab208117) ab208117 staining Caldesmon/CDM in HeLa cells. The cells were fixed with 4% formaldehyde (10 min), permeabilized with 0.1% Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at +4°C with ab208117 at 1/100 dilution (shown in red) and **ab195887**, Mouse monoclonal to alpha Tubulin (Alexa Fluor<sup>®</sup> 488), at 1/250 dilution (shown in green). Nuclear DNA was labelled with DAPI (shown in blue).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

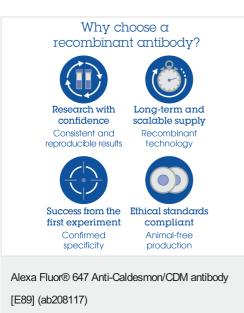

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

#### Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise,

please visit https://www.abcam.co.jp/abpromise or contact our technical team.

#### **Terms and conditions**

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors