abcam

Product datasheet

Alexa Fluor® 647 Anti-CNPase antibody [11-5B] ab201679

画像数 1

製品の概要

特記事項

製品名 Alexa Fluor® 647 Anti-CNPase antibody [11-5B]

製品の詳細 Alexa Fluor® 647 Mouse monoclonal [11-5B] to CNPase

由来種 Mouse

標識 Alexa Fluor® 647. Ex: 652nm, Em: 668nm

アプリケーション 適用あり: IHC-P 種交差性 交差種: Human

交差が予測される動物種: Mouse, Rat, Sheep, Rabbit, Cow, Dog, Pig, Rhesus monkey 4 非

交差種: Chicken, Guinea pig

免疫原 Full length native protein (purified) corresponding to Human CNPase.

ポジティブ・コントロール IHC-FFPE sections - Human cerebral cortex (normal)

> Alexa Fluor® is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor® dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor® dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor[®] dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@thermofisher.com.

> The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

バッファー pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: PBS, 30% Glycerol (glycerin, glycerine), 1% BSA

Batches contain 0.4M arginine.

精製度 Affinity purified

ポリ/モノ モノクローナル

ウローン名 11-5B **アイソタイプ I**gG1

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab201679の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		1/100. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.


ターゲット情報

配列類似性 Belongs to the cyclic nucleotide phosphodiesterase family.

細胞内局在 Membrane. Melanosome. Firmly bound to membrane structures of brain white matter. Identified

by mass spectrometry in melanosome fractions from stage I to stage ${\sf IV}$.

画像

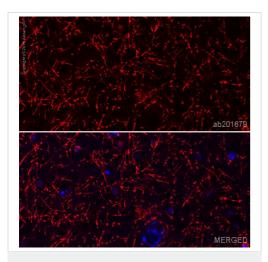

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Alexa Fluor® 647 Anti-CNPase antibody [11-5B] (ab201679)

IHC image of CNPase staining in a section of formalin-fixed paraffin-embedded normal human cerebral cortex.

The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6) in a Dako Pascal pressure cooker using the standard factory-set regime. Non-specific protein-protein interactions were then blocked in TBS containing 0.025% (v/v) Triton X-100, 0.3M (w/v) glycine and 1% (w/v) BSA for 1h at room temperature. The section was then incubated overnight at +4°C in TBS containing 0.025% (v/v) Triton X-100 and 1% (w/v) BSA with ab201679 at 1/100 (shown in red). Nuclear DNA was labelled with DAPI (shown in blue). The section was then mounted using Fluoromount[®].

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

For other IHC staining systems (automated and non-automated), customers should optimize variable parameters such as antigen retrieval conditions, antibody concentrations and incubation times.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors