abcam

Product datasheet

Alexa Fluor® 647 Anti-Ki67 antibody [EPR3610] ab196907

КО 評価済 มาวชาวง RabMAb

★★★★★ <u>1 Abreviews</u> <u>2 References</u> 画像数 4

製品の概要

製品名	Alexa Fluor® 647 Anti-Ki67 antibody [EPR3610]		
製品の詳細	Alexa Fluor® 647 Rabbit monoclonal [EPR3610] to Ki67		
由来種	Rabbit		
標識	Alexa Fluor® 647. Ex: 652nm, Em: 668nm		
アプリケーション	適用あり: ICC/IF		
種交差性	交差種: Human		
免疫原	Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.		
ポジティブ・コントロール	ICC: HeLa and wild-type HAP1 cells.		
特記事項	 This product is a recombinant monoclonal antibody, which offers several advantages including: High batch-to-batch consistency and reproducibility Improved sensitivity and specificity Long-term security of supply Animal-free production For more information see here. Our RabMAb[®] technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to RabMAb[®] patents. Alexa Fluor[®] is a registered trademark of Molecular Probes, Inc, a Thermo Fisher Scientific Company. The Alexa Fluor[®] dye included in this product is provided under an intellectual property license from Life Technologies Corporation. As this product contains the Alexa Fluor[®] dye, the purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). As this product contains the Alexa Fluor[®] dye the sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: in manufacturing; (ii) to provide a service, information, or data in return for payment (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are sold for use in research. For information on purchasing a license to this product for purposes other than research, contact Life Technologies Corporation, 5781 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@thermofisher.com. 		

製品の特性	
製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C. Avoid freeze / thaw cycle. Store In the Dark.
バッファー	pH: 7.40 Preservative: 0.02% Sodium azide Constituents: PBS, 30% Glycerol (glycerin, glycerine), 1% BSA
精製度	Protein A purified
ポリ/モノ	モノクローナル
クローン名	EPR3610
アイソタイプ	lgG

アプリケーション

The Abpromise guaranteeAbpromise保証は、
次のテスト済みアプリケーションにおけるab196907の使用に適用されますアプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
ICC/IF	★ ★ ★ ★ ★ (1)	1/100. If fixing cells in 4% PFA, it is recommended to permeabilized cells with 0.1% Triton-X for 5 min.

ターゲット情報


機能	Required to maintain individual mitotic chromosomes dispersed in the cytoplasm following nuclear envelope disassembly (PubMed:27362226). Associates with the surface of the mitotic chromosome, the perichromosomal layer, and covers a substantial fraction of the chromosome surface (PubMed:27362226). Prevents chromosomes from collapsing into a single chromatin mass by forming a steric and electrostatic charge barrier: the protein has a high net electrical charge and acts as a surfactant, dispersing chromosomes and enabling independent chromosome motility (PubMed:27362226). Binds DNA, with a preference for supercoiled DNA and AT-rich DNA (PubMed:10878551). Does not contribute to the internal structure of mitotic chromosomes (By similarity). May play a role in chromatin organization (PubMed:24867636). It is however unclear whether it plays a direct role in chromatin organization or whether it is an indirect consequence of its function in maintaining mitotic chromosomes dispersed.
配列類似性	Contains 1 FHA domain. Contains 16 K167R repeats. Contains 1 PP1-binding domain.
発生段階	Expression occurs preferentially during late G1, S, G2 and M phases of the cell cycle, while in cells in G0 phase the antigen cannot be detected (at protein level) (PubMed:6206131). Present at highest level in G2 phase and during mitosis (at protein level). In interphase, forms fiber-like structures in fibrillarin-deficient regions surrounding nucleoli (PubMed:2674163, PubMed:8799815).
翻訳後修飾	Phosphorylated. Hyperphosphorylated in mitosis (PubMed:10502411, PubMed:10653604).

細胞内局在

Hyperphosphorylated form does not bind DNA.

Chromosome. Nucleus. Nucleus, nucleolus. Associates with the surface of the mitotic chromosome, the perichromosomal layer, and covers a substantial fraction of the mitotic chromosome surface (PubMed:27362226). Associates with satellite DNA in G1 phase (PubMed:9510506). Binds tightly to chromatin in interphase, chromatin-binding decreases in mitosis when it associates with the surface of the condensed chromosomes (PubMed:15896774, PubMed:22002106). Predominantly localized in the G1 phase in the perinucleolar region, in the later phases it is also detected throughout the nuclear interior, being predominantly localized in the nuclear matrix (PubMed:22002106).

画像

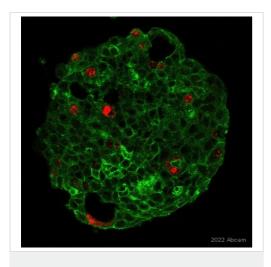

Immunocytochemistry - Alexa Fluor® 647 Anti-Ki67 antibody [EPR3610] (ab196907)


This image is courtesy of an Abreview submitted by Tania Cruz

Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 647 Anti-Ki67 antibody [EPR3610] (ab196907)

Immunocytochemistry analysis of buffered formalin-fixed 0.05% Tween 20 permeabilized human rectal cancer spheroids, SW837 staining with ab196907 at 1/100 dilution. Samples were incubated with the primary antibody with Blocking solution for 12 hours at 4°C. Blocking was done using 1% BSA + 10% FBS for 1 hour at 25°C.

ab196907 staining Ki67 in wild-type HAP1 cells (top panel) and Ki67 knockout HAP1 cells (bottom panel). The cells were fixed with 100% methanol (5min), permeabilized with 0.1% Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated with ab196907 at 1/100 dilution (shown in red) and ab195887 at 1/250 dilution (shown in green) overnight at +4°C. Nuclear DNA was labelled in blue with DAPI.

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

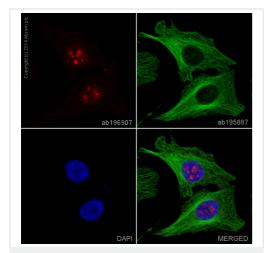

Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 647 Anti-Ki67 antibody [EPR3610] (ab196907)



(ab196907)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise,

ab196907 staining Ki67 in HeLa cells. The cells were fixed with 4% formaldehyde (10 min), permeabilized in 0.1% Triton X-100 for 5 minutes and then blocked in 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated with ab196907 at 1/100 dilution (shown in red) and **ab195887**, Mouse monoclonal [DM1A] to alpha Tubulin (Alexa Fluor[®] 488, shown in green) at 1/167 dilution, overnight at +4°C. Nuclear DNA was labelled in blue with DAPI.

This product also gave a positive signal under the same testing conditions in HeLa cells fixed with 100% methanol (5 min).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8)

please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors