

### Anti-ATRX antibody [AMab-6] ab264537

リコンビナント

画像数 2

#### 製品の概要

製品名	Anti-ATRX antibody [AMab-6]
製品の詳細	Mouse monoclonal [AMab-6] to ATRX
由来種	Mouse
アプリケーション	<b>適用あり:</b> IHC-P <b>適用なし:</b> ICC,IHC-Fr or WB
種交差性	<b>交差種:</b> Human
免疫原	Recombinant fragment. This information is proprietary to Abcam and/or its suppliers.
ポジティブ・コントロール	IHC-P: IDH-mutant human oligodendroglioma, IDH-mutant human diffuse astrocytoma, Human breast carcinoma.
特記事項	<p>This antibody clone is manufactured by Abcam. If you require a custom buffer formulation or conjugation for your experiments, please contact <a href="mailto:orders@abcam.com">orders@abcam.com</a>.</p> <p>This product is a recombinant monoclonal antibody, which offers several advantages including:</p> <ul style="list-style-type: none"><li>- High batch-to-batch consistency and reproducibility</li><li>- Improved sensitivity and specificity</li><li>- Long-term security of supply</li><li>- Animal-free production</li></ul> <p>For more information <a href="#">see here</a>.</p>

#### 製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C. Avoid freeze / thaw cycle.
バッファー	pH: 7.20 Preservative: 0.01% Sodium azide Constituents: 59% PBS, 40% Glycerol (glycerin, glycerine), 0.05% BSA
精製度	Protein A purified
ポリ/モノ	モノクローナル
クローン名	AMab-6
アイソタイプ	IgG1

## アプリケーション

**The Abpromise guarantee**      **Abpromise保証は、** 次のテスト済みアプリケーションにおけるab264537の使用に適用されます  
アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		1/10000. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

**追加情報**      Is unsuitable for ICC,IHC-Fr or WB.

## ターゲット情報

**機能**      Involved in transcriptional regulation and chromatin remodeling. Facilitates DNA replication in multiple cellular environments and is required for efficient replication of a subset of genomic loci. Binds to DNA tandem repeat sequences in both telomeres and euchromatin and in vitro binds DNA quadruplex structures. May help stabilizing G-rich regions into regular chromatin structures by remodeling G4 DNA and incorporating H3.3-containing nucleosomes. Catalytic component of the chromatin remodeling complex ATRX:DAXX which has ATP-dependent DNA translocase activity and catalyzes the replication-independent deposition of histone H3.3 in pericentric DNA repeats outside S-phase and telomeres, and the in vitro remodeling of H3.3-containing nucleosomes. Its heterochromatin targeting is proposed to involve a combinatorial readout of histone H3 modifications (specifically methylation states of H3K9 and H3K4) and association with CBX5. Involved in maintaining telomere structural integrity in embryonic stem cells which probably implies recruitment of CBX5 to telomers. Reports on the involvement in transcriptional regulation of telomeric repeat-containing RNA (TERRA) are conflicting; according to a report, it is not sufficient to decrease chromatin condensation at telomers nor to increase expression of telomeric RNA in fibroblasts (PubMed:24500201). May be involved in telomere maintenance via recombination in ALT (alternative lengthening of telomeres) cell lines. Acts as negative regulator of chromatin incorporation of transcriptionally repressive histone H2AFY, particularly at telomeres and the alpha-globin cluster in erythroleukemic cells. Participates in the allele-specific gene expression at the imprinted IGF2/H19 gene locus. On the maternal allele, required for the chromatin occupancy of SMC1 and CTCF within the H19 imprinting control region (ICR) and involved in establishment of histone tails modifications in the ICR. May be involved in brain development and facial morphogenesis. Binds to zinc-finger coding genes with atypical chromatin signatures and regulates its H3K9me3 levels. Forms a complex with ZNF274, TRIM28 and SETDB1 to facilitate the deposition and maintenance of H3K9me3 at the 3' exons of zinc-finger genes (PubMed:27029610).

**組織特異性**      Ubiquitous.

**関連疾患**      Alpha-thalassemia mental retardation syndrome, X-linked  
Mental retardation, X-linked, syndromic, with hypotonic facies 1  
Alpha-thalassemia myelodysplasia syndrome

**配列類似性**      Belongs to the SNF2/RAD54 helicase family.  
Contains 1 ADD domain.  
Contains 1 GATA-type zinc finger.

Contains 1 helicase ATP-binding domain.

Contains 1 helicase C-terminal domain.

Contains 1 PHD-type zinc finger.

## ドメイン

The ADD domain predominantly interacts with histone H3 trimethylated at 'Lys-10' (H3K9me3) (and to a lesser extent H3 mono-or dimethylated at 'Lys-10') and simultaneously to histone H3 unmethylated at 'Lys-5' (H3K4me0). The interaction with H3K9me3 is disrupted by the presence of H3K4me3 suggesting a readout of the combined histone H3 methylation state.

Contains one Pro-Xaa-Val-Xaa-Leu (PxVxL) motif, which is required for interaction with chromoshadow domains. This motif requires additional residues -7, -6, +4 and +5 of the central Val which contact the chromoshadow domain.


## 翻訳後修飾

Phosphorylated at serine residues during mitose. Phosphorylation may promote the release from the nuclear matrix and progression to mitosis.

## 細胞内局在

Nucleus. Chromosome, telomere. Nucleus, PML body. Associated with pericentromeric heterochromatin during interphase and mitosis, probably by interacting with CBX5/HP1 alpha. Colocalizes with histone H3.3, DAXX, HIRA and ASF1A at PML-nuclear bodies. Colocalizes with cohesin (SMC1 and SMC3) and MECP2 at the maternal H19 ICR (By similarity).

## 画像

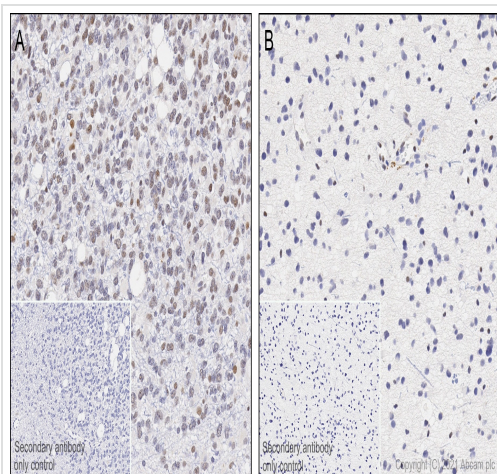


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ATRX antibody [AMab-6] (ab264537)

Immunohistochemical analysis of paraffin-embedded Image A IDH-mutant human oligodendroglioma; Image B IDH-mutant human diffuse astrocytoma tissue labeling None with abab264537 at 1/10000 dilution followed by a ready to use LeicaDS9800 (Bond™ Polymer Refine Detection) was used. Image A: nuclear staining on tumor cells without ATRX mutation in IDH-mutant human oligodendroglioma; Image B: negative staining of ATRX in the nucleus of tumor cells with ATRX mutation, whereas positive staining in the vascular endothelial cells without ATRX mutation (PMID: 27788029). The section was incubated with ab264537 for 15 mins at room temperature. The immunostaining was performed on a Leica Biosystems BOND® RX instrument Counterstained with Hematoxylin.

Secondary antibody only control: Secondary antibody is a ready to use LeicaDS9800 (Bond™ Polymer Refine Detection) was used.

Heat mediated antigen retrieval with Citrate buffer (pH 6.0, epitope retrieval solution 1) for 20 mins

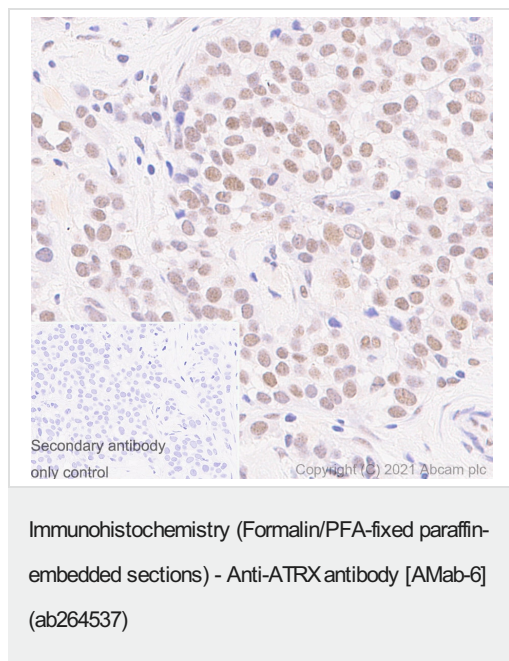

Immunohistochemical analysis of paraffin-embedded Human breast carcinoma tissue labeling None with abab264537 at 1/10000 dilution followed by a ready to use LeicaDS9800 (Bond™ Polymer Refine Detection) was used. Nuclear staining on human breast carcinoma. The section was incubated with ab264537 for 15 mins at room temperature. The immunostaining was performed on a Leica Biosystems BOND® RX instrument Counterstained with Hematoxylin.

Secondary antibody only control: Secondary antibody is a ready to use LeicaDS9800 (Bond™ Polymer Refine Detection) was used.

Heat mediated antigen retrieval with Citrate buffer (pH 6.0, epitope retrieval solution 1) for 20 mins

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

### Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

### Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors