

Anti-CD45 antibody [F10-89-4] - Hematopoietic Stem Cell Marker ab30470

★★★★★ [1 Abreviews](#) [10 References](#) [画像数 4](#)

製品の概要

製品名	Anti-CD45 antibody [F10-89-4] - Hematopoietic Stem Cell Marker
製品の詳細	Mouse monoclonal [F10-89-4] to CD45 - Hematopoietic Stem Cell Marker
由来種	Mouse
特異性	Clone F10-89-4 reacts with all forms of CD45 expressed by all haematopoietic cells, except erythrocytes, having a higher level of expression on lymphocytes than on granulocytes.
アプリケーション	適用あり: ICC/IF, Flow Cyt, IHC-P
種交差性	交差種: Human
免疫原	Tissue, cells or virus corresponding to Human CD45. Human T lymphocytes.
ポジティブ・コントロール	IHC-P: Human tonsil (normal) tissue. ICC/IF: Jurkat cells. Flow Cyt: Jurkat cells, whole blood.
特記事項	<p>This antibody clone is manufactured by Abcam. If you require a custom buffer formulation or conjugation for your experiments, please contact orders@abcam.com.</p> <p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As</p>

製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long term. Avoid freeze / thaw cycle.
バッファー	<p>pH: 7.40</p> <p>Preservative: 0.02% Sodium azide</p> <p>Constituent: PBS</p>
精製度	Protein G purified
一次抗体 備考	Clone F10-89-4 reacts with all forms of CD45 expressed by all haematopoietic cells, except

	erythrocytes, having a higher level of expression on lymphocytes than on granulocytes.
ポリ/モノ	モノクローナル
クローン名	F10-89-4
ミエローマ	NS1
アイソタイプ	IgG2a


アプリケーション

The Abpromise guarantee Abpromise保証は、 次のテスト済みアプリケーションにおけるab30470の使用に適用されます
アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
ICC/IF	★★★★★ (1)	Use a concentration of 5 - 10 µg/ml.
Flow Cyt		1/10 - 1/50. ab170191 - Mouse monoclonal IgG2a, is suitable for use as an isotype control with this antibody.
IHC-P		Use a concentration of 1 - 5 µg/ml. Perform heat mediated antigen retrieval before commencing with IHC staining protocol.

ターゲット情報

機能	Protein tyrosine-protein phosphatase required for T-cell activation through the antigen receptor. Acts as a positive regulator of T-cell coactivation upon binding to DPP4. The first PTPase domain has enzymatic activity, while the second one seems to affect the substrate specificity of the first one. Upon T-cell activation, recruits and dephosphorylates SKAP1 and FYN.
関連疾患	Defects in PTPRC are a cause of severe combined immunodeficiency autosomal recessive T-cell-negative/B-cell-positive/NK-cell-positive (T(-)B(+)NK(+) SCID) [MIM:608971]. A form of severe combined immunodeficiency (SCID), a genetically and clinically heterogeneous group of rare congenital disorders characterized by impairment of both humoral and cell-mediated immunity, leukopenia, and low or absent antibody levels. Patients present in infancy recurrent, persistent infections by opportunistic organisms. The common characteristic of all types of SCID is absence of T-cell-mediated cellular immunity due to a defect in T-cell development. Genetic variations in PTPRC are involved in multiple sclerosis susceptibility (MS) [MIM:126200]. MS is a neurodegenerative disorder characterized by the gradual accumulation of focal plaques of demyelination particularly in the periventricular areas of the brain. Peripheral nerves are not affected. Onset usually in third or fourth decade with intermittent progression over an extended period. The cause is still uncertain.
配列類似性	Belongs to the protein-tyrosine phosphatase family. Receptor class 1/6 subfamily. Contains 2 fibronectin type-III domains. Contains 2 tyrosine-protein phosphatase domains.
ドメイン	The first PTPase domain interacts with SKAP1.
翻訳後修飾	Heavily N- and O-glycosylated.
細胞内局在	Membrane. Membrane raft. Colocalized with DPP4 in membrane rafts.

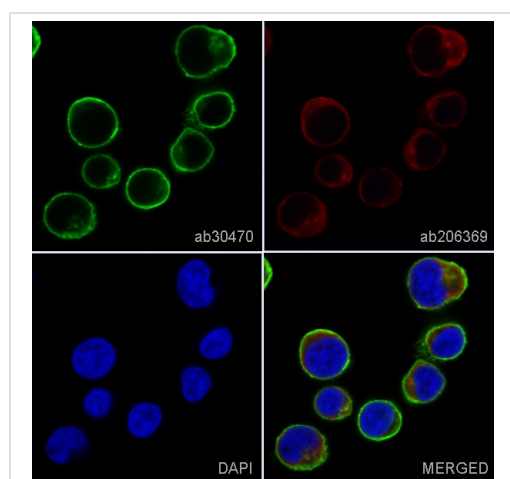


Immunocytochemistry/ Immunofluorescence - Anti-CD45 antibody [F10-89-4] - Hematopoietic Stem Cell Marker (ab30470)

ab30470 staining CD45 in Jurkat (Human T cell leukemia cell line from peripheral blood) cells. The cells were fixed with 4% formaldehyde (10 minutes), then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1 hour. The cells were then incubated overnight at +4°C with ab30470 at 5 µg/ml (shown in green) and **ab206369**, Rabbit monoclonal to beta Tubulin (Alexa Fluor® 594), at 1/250 dilution (shown in red). This was followed by an incubation at room temperature for 1 hour with **ab150117**, Goat Anti-Mouse IgG H&L (Alexa Fluor® 488) preadsorbed, at 1 µg/ml (shown in green). Nuclear DNA was labelled with DAPI (shown in blue).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

This product also gave a positive signal under the same testing conditions in Jurkat cells fixed with 80% methanol (5 minutes).

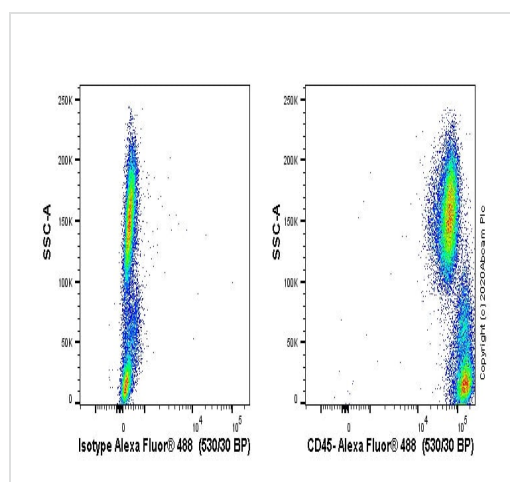


Flow Cytometry - Anti-CD45 antibody [F10-89-4] - Hematopoietic Stem Cell Marker (ab30470)

Flow cytometry staining of human whole blood with ab30470 (right) or mouse IgG2ak; (**ab18413**) isotype (left). Red blood cells of 200 µL blood were lysed, then cells were incubated for 30 min on ice in 1x PBS containing 10 µg/ml human IgG and 10 % normal goat serum to block FC receptors and non-specific protein-protein interaction followed by the antibody (ab30470) or mouse IgG2ak; (**ab18413**) isotype (1×10^6 in 100 µL at 0.2 µg/ml) for 30 min on ice.


The secondary antibody Goat anti-mouse IgG H&L (Alexa Fluor® 488, pre-adsorbed) (**ab150117**) was used at 1/2000 dilution for 30 min on ice.

Acquisition of >30000 events were collected using a 50 mW Blue laser (488nm) and 530/30 bandpass filter. Events were gated on alive cells.

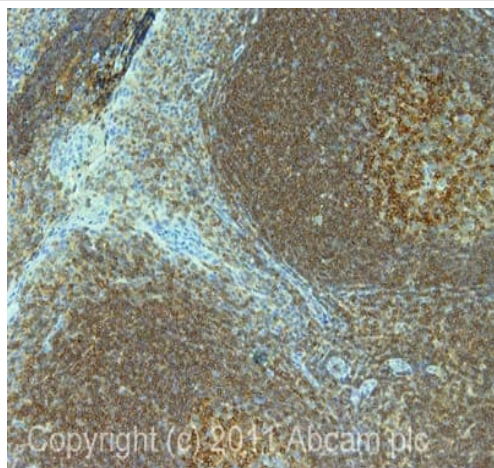

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-CD45 antibody [F10-89-4]
- Hematopoietic Stem Cell Marker (ab30470)

IHC image of CD45 staining in human tonsil formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6, epitope retrieval solution 1) for 20 minutes. The section was then incubated with ab30470, 1 µg/ml, for 15 minutes at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

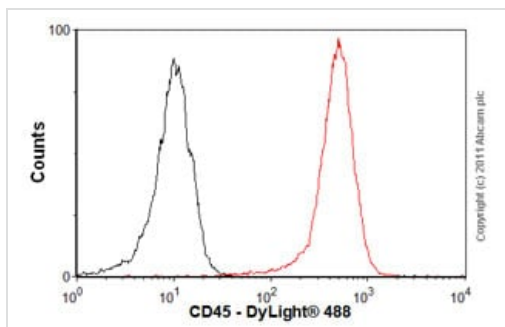

Flow Cytometry - Anti-CD45 antibody [F10-89-4] - Hematopoietic Stem Cell Marker (ab30470)

Overlay histogram showing Jurkat (Human T cell leukemia cell line from peripheral blood) cells stained with ab30470 (red line). The cells were fixed with 4% paraformaldehyde (10 minutes) and incubated in 1x PBS / 10% normal goat serum / 0.3M glycine to block non-specific protein-protein interactions. The cells were then incubated with the antibody (ab30470, 1 µg/1 x 10⁶ cells) for 30 minutes at 22°C. The secondary antibody used was DyLight® 488 goat anti-mouse IgG (H+L) ([ab96879](#)) at 1/500 dilution for 30 minutes at 22°C. Isotype control antibody (black line) was mouse IgG2a [ICIGG2A] ([ab91361](#), 1 µg/1 x 10⁶ cells) used under the same conditions. Acquisition of >5,000 events was performed. This antibody gave a positive signal in Jurkat cells fixed with 80% methanol (5 minutes) used under the same conditions.

Please note that Abcam do not have any data for use of this antibody on non-fixed cells. We welcome any customer feedback.

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours

- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors