# abcam

## Product datasheet

## Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free ab269369



יעדער RabMAb

## 画像数 5

## 製品の概要

製品名 Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free

製品の詳細 Rabbit monoclonal [EPR23222-24] to Eph receptor B4/HTK - BSA and Azide free

由来種 Rabbit

アプリケーション 適用あり: IHC-P, WB

適用なし: Flow Cyt,ICC/IF or IP

種交差性 交差種: Mouse, Human

免疫原 Recombinant fragment. This information is proprietary to Abcam and/or its suppliers.

ポジティブ・コントロール WB: HEK293T, HT-29, HUVEC, HCT116, MCF7, T47D, PC-3, NIH:OVCAR-3, SK-BR-3, GR-M,

4T1 and NIH/3T3 lysates. IHC-P: Human colon, Human colon cancer and Human breast cancer

tissues.

特記事項 ab269369 is the carrier-free version of ab254301.

> Our carrier-free antibodies are typically supplied in a PBS-only formulation, purified and free of BSA, sodium azide and glycerol. The carrier-free buffer and high concentration allow for increased conjugation efficiency.

This conjugation-ready format is designed for use with fluorochromes, metal isotopes, oligonucleotides, and enzymes, which makes them ideal for antibody labelling, functional and cellbased assays, flow-based assays (e.g. mass cytometry) and Multiplex Imaging applications.

Use our conjugation kits for antibody conjugates that are ready-to-use in as little as 20 minutes with <1 minute hands-on-time and 100% antibody recovery: available for fluorescent dyes, HRP, biotin and gold.

This product is compatible with the Maxpar® Antibody Labeling Kit from Fluidigm, without the need for antibody preparation. Maxpar<sup>®</sup> is a trademark of Fluidigm Canada Inc.

This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility
- Improved sensitivity and specificity
- Long-term security of supply
- Animal-free production

For more information see here.

Our RabMAb® technology is a patented hybridoma-based technology for making rabbit

#### 製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C. Do Not Freeze.

**バッファー** pH: 7.2

Constituent: PBS

キャリア・フリー はい

精製度 Protein A purified

**ポリ**(モノ モノクローナル **ウローン名** EPR23222-24

アイソタイプ IgG

#### アプリケーション

**The Abpromise guarantee** <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab269369の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P		Use at an assay dependent concentration. Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.
WB		Use at an assay dependent concentration. Detects a band of approximately 120 kDa (predicted molecular weight: 108 kDa).

追加情報 Is unsuitable for Flow Cyt,ICC/IF or IP.

Autophosphorylated.

Membrane.

## ターゲット情報


翻訳後修飾

細胞内局在

機能
Receptor for members of the ephrin-B family. Binds to ephrin-B2. May have a role in events mediating differentiation and development.

Abundantly expressed in placenta and in a range of primary tissues and malignant cell lines. Expressed in fetal, but not adult, brain, and in primitive and myeloid, but not lymphoid, hematopoietic cells.

配列類似性
Belongs to the protein kinase superfamily. Tyr protein kinase family. Ephrin receptor subfamily. Contains 2 fibronectin type-III domains.
Contains 1 protein kinase domain.
Contains 1 SAM (sterile alpha motif) domain.

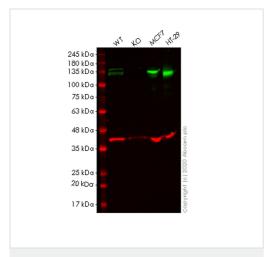

Western blot - Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free (ab269369)

**All lanes :** Anti-Eph receptor B4/HTK antibody [EPR23222-24] (ab254301) at 1/1000 dilution

Lane 1: Wild-type HEK293T cell lysate

Lane 2: EPHB4 knockout HEK293T cell lysate

Lane 3 : MCF7 cell lysate

Lane 4 : HT-29 cell lysate

Lysates/proteins at 20 µg per lane.

## Secondary


**All lanes :** Goat anti-Rabbit lgG H&L (IRDye® 800CW) preadsorbed (ab216773) at 1/10000 dilution

Predicted band size: 108 kDa
Observed band size: 125 kDa

This data was developed using the same antibody clone in a different buffer formulation (<u>ab254301</u>).

**Lanes 1-4:** Merged signal (red and green). Green - <u>ab254301</u> observed at 125 kDa. Red - loading control <u>ab8245</u> observed at 36 kDa.

ab254301 Anti-Eph receptor B4/HTK antibody [EPR23222-24] was shown to specifically react with Eph receptor B4/HTK in wild-type HEK293T cells. Loss of signal was observed when knockout cell line ab266733 (knockout cell lysate ab257429) was used. Wild-type and Eph receptor B4/HTK knockout samples were subjected to SDS-PAGE. ab254301 and Anti-GAPDH antibody [6C5] - Loading Control (ab8245) were incubated overnight at 4°C at 1 in 1000 dilution and 1 in 20000 dilution respectively. Blots were developed with Goat anti-Rabbit IgG H&L (IRDye® 800CW) preadsorbed (ab216773) and Goat anti-Mouse IgG H&L (IRDye® 680RD) preadsorbed (ab216776) secondary antibodies at 1 in 20000 dilution for 1 hour at room temperature before imaging.

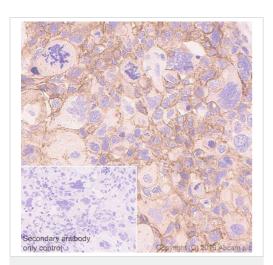


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free (ab269369)

Immunohistochemical analysis of paraffin-embedded Human breast cancer tissue labeling Eph receptor B4 with <u>ab254301</u> at 1/1000 dilution (0.545ug/ml) followed by a Rabbit specific IHC polymer detection kit HRP/DAB (<u>ab209101</u>) at a ready to use dilution. Mainly membranous staining in cancer cells of human breast cancer (PMID: 16816380) is observed. The section was incubated with <u>ab255611</u> for 30 mins at room temperature. The immunostaining was performed on a Leica Biosystems BOND<sup>®</sup> RX instrument. Counterstained with Hematoxylin.

Secondary antibody only control: Secondary antibody is Rabbit specific IHC polymer detection kit HRP/DAB (<u>ab209101</u>) at Ready to use dilution.

Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0, epitope retrieval solution 2) for 20 mins.

This data was developed using the same antibody clone in a different buffer formulation containing PBS, BSA, glycerol, and sodium azide (ab254301).

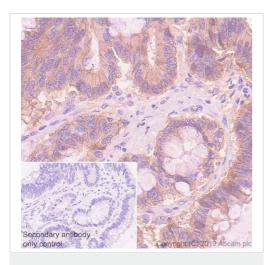

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free (ab269369)

Immunohistochemical analysis of paraffin-embedded Human colon cancer tissue labeling Eph receptor B4 with <a href="mailto:ab254301">ab254301</a> at 1/1000 dilution (0.545ug/ml) followed by a Rabbit specific IHC polymer detection kit HRP/DAB (<a href="mailto:ab209101">ab209101</a>) at a ready to use dilution. Positive staining in cancer cells of human colon cancer (PMID: 19366806) is observed. The section was incubated with <a href="mailto:ab255611">ab255611</a> for 30 mins at room temperature. The immunostaining was performed on a Leica Biosystems BOND® RX instrument. Counterstained with Hematoxylin.

Secondary antibody only control: Secondary antibody is Rabbit specific IHC polymer detection kit HRP/DAB (ab209101) at Ready to use dilution.

Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0, epitope retrieval solution 2) for 20 mins.

This data was developed using the same antibody clone in a different buffer formulation containing PBS, BSA, glycerol, and sodium azide (ab254301).



Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Eph receptor B4/HTK antibody [EPR23222-24] - BSA and Azide free (ab269369)

Immunohistochemical analysis of paraffin-embedded Human colon tissue labeling Eph receptor B4 with <u>ab254301</u> at 1/1000 dilution (0.545ug/ml) followed by a Rabbit specific IHC polymer detection kit HRP/DAB (<u>ab209101</u>) at a ready to use dilution. Weak positive staining in base of crypts of human colon (PMID:19366806) is observed. The section was incubated with <u>ab255611</u> for 30 mins at room temperature. The immunostaining was performed on a Leica Biosystems BOND® RX instrument. Counterstained with Hematoxylin.

Secondary antibody only control: Secondary antibody is Rabbit specific IHC polymer detection kit HRP/DAB (<u>ab209101</u>) at Ready to use dilution.

Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0, epitope retrieval solution 2) for 20 mins.

This data was developed using the same antibody clone in a different buffer formulation containing PBS, BSA, glycerol, and sodium azide (ab254301).



BSA and Azide free (ab269369)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

## Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish

- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.co.jp/abpromise">https://www.abcam.co.jp/abpromise</a> or contact our technical team.

#### Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors