

Anti-Giantin antibody - Golgi Marker ab80864

KO 評価済

★★★★☆ 2 Abreviews 24 References 画像数 4

製品の概要

製品名	Anti-Giantin antibody - Golgi Marker
製品の詳細	Rabbit polyclonal to Giantin - Golgi Marker
由来種	Rabbit
アプリケーション	適用あり: IHC-P, ICC/IF, WB
種交差性	交差種: Human
免疫原	Synthetic peptide corresponding to Human Giantin aa 1250-1350 conjugated to keyhole limpet haemocyanin. (Peptide available as ab101412)
ポジティブ・コントロール	ICC/IF: MCF7 cells. IHC-P: Human kidney tissue.
特記事項	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As</p>

製品の特性

製品の状態	Liquid
保存方法	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -80°C. Avoid freeze / thaw cycle.
バッファー	<p>pH: 7.40</p> <p>Preservative: 0.02% Sodium azide</p> <p>Constituent: PBS</p> <p>Batches of this product that have a concentration < 1mg/ml may have BSA added as a stabilising agent. If you would like information about the formulation of a specific lot, please contact our scientific support team who will be happy to help.</p>
精製度	Immunogen affinity purified

ポリクローナル
IgG

The Abpromise guarantee Abpromise保証は、次のテスト済みアプリケーションにおけるab80864の使用に適用されます
アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P	★★★★☆ (1)	Use a concentration of 0.1 - 1 µg/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.
ICC/IF	★★★★☆ (1)	Use a concentration of 1 µg/ml.
WB		1/1000. Predicted molecular weight: 376 kDa.


機能	May participate in forming intercisternal cross-bridges of the Golgi complex.
細胞内局在	Golgi apparatus membrane.

Immunofluorescence images of cells stained for Ab7291 (red), Ab80864 (green), and DAPI (blue). The top row shows Ab7291 (left) and Ab80864 (right). The bottom row shows DAPI (left) and the merged image (right).

ab80864 staining Giantin in HeLa cells. The cells were fixed with 100% methanol (5 min), permeabilized with 0.1% PBS-Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at 4°C with ab80864 at 1µg/ml and **ab7291**, Mouse monoclonal [DM1A] to alpha Tubulin - Loading Control. Cells were then incubated with **ab150081**, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 488), pre-adsorbed at 1/1000 dilution (shown in green) and **ab150120**, Goat polyclonal Secondary Antibody to Mouse IgG - H&L (Alexa Fluor® 594), pre-adsorbed at 1/1000 dilution (shown in pseudocolour red). Nuclear DNA was labelled with DAPI (shown in blue).

Also suitable in cells fixed with 4% paraformaldehyde (10 min).

Image was acquired with a high-content analyser (Operetta CLS, Perkin Elmer) and a maximum intensity projection of confocal sections is shown.

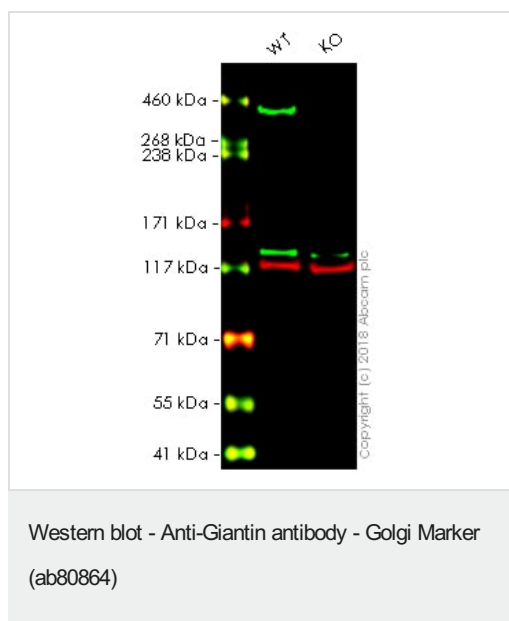

All lanes : Anti-Giantin antibody - Golgi Marker (ab80864) at 1/1000 dilution

Lane 1 : Wild-type HAP1 whole cell lysate


Lane 2 : GOLGB1 knockout HAP1 whole cell lysate

Lysates/proteins at 20 µg per lane.

Predicted band size: 376 kDa

Lanes 1 - 2: Merged signal (red and green). Green - ab80864 observed at 377kDa. Red - loading control, [ab130007](#), observed at 130kDa.

ab80864 was shown to recognize GOLGB1 in wild-type HAP1 cells as signal was lost at the expected MW in GOLGB1 knockout cells. Additional cross-reactive bands were observed in the wild-type and knockout cells. Wild-type and GOLGB1 knockout samples were subjected to SDS-PAGE. Ab80864 and [ab130007](#) (Mouse anti-Vinculin loading control) were incubated overnight at 4°C at 1/1000 dilution and 1/20000 dilution respectively. Blots were developed with Goat anti-Rabbit IgG H&L (IRDye® 800CW) preabsorbed [ab216773](#) and Goat anti-Mouse IgG H&L (IRDye® 680RD) preabsorbed [ab216776](#) secondary antibodies at 1/20000 dilution for 1 hour at room temperature before imaging.

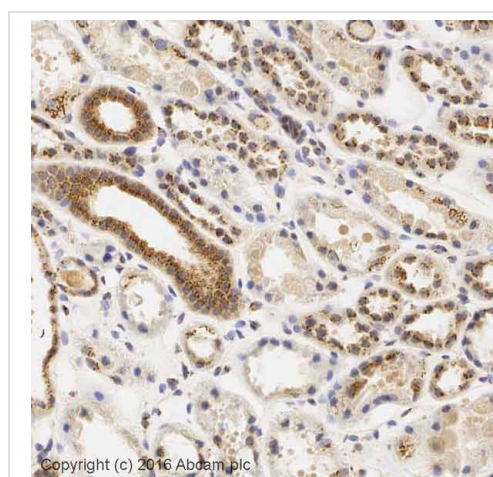


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Giantin antibody - Golgi Marker (ab80864)

IHC image of ab80864 staining in human kidney formalin fixed paraffin embedded tissue section*, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab80864, 0.1µg/ml for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

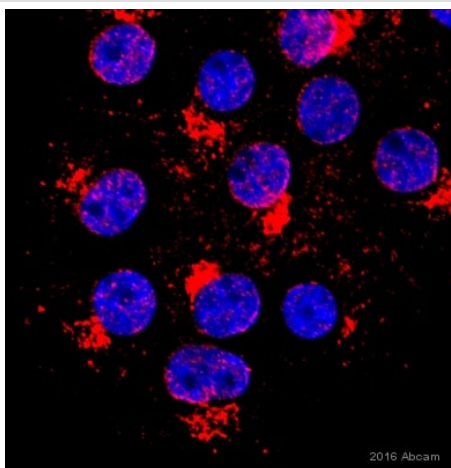

Ab80864 staining Giantin in WIF-B cells (hepatoma-derived hybrid cell line) by ICC/IF (Immunocytochemistry/Immunofluorescence). Cells were fixed with formaldehyde; permeabilized with 0.2% Triton X-100 and blocked with 1% Donkey serum in 0.1% PBST for 1 hour at 22°C. Samples were incubated at 1/25 dilution for 3 hours at 22°C. An Alexa Fluor® 594 Donkey anti rabbit was used as the secondary antibody at 1/200 dilution.

Immunocytochemistry/ Immunofluorescence - Anti-Giantin antibody - Golgi Marker (ab80864)

This image is courtesy of an anonymous Abreview

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors