abcam

Product datasheet

Anti-Histone H2A (acetyl K5) antibody ab1764

★★★★★ 10 Abreviews 12 References 画像数 1

製品の概要

製品名 Anti-Histone H2A (acetyl K5) antibody

製品の詳細 Rabbit polyclonal to Histone H2A (acetyl K5)

由来種 Rabbit

 アプリケーション
 適用あり: IHC-P

 種交差性
 交差種: Human

免疫原 Synthetic peptide corresponding to Histone H2A (acetyl K5) conjugated to ovalbumin.

■ Run BLAST with EXPASY ■ Run BLAST with S NCBI

特記事項 <u>ab171870</u> - Rabbit polyclonal IgG, is suitable for use as an isotype control with this antibody.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

パッファー Preservative: 0.02% Sodium azide

精製度 Whole antiserum ポリ/モノ ポリクローナル

アイソタイプ lqG

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab1764の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

1

アプリケーション	Abreviews	特記事項
IHC-P	**** <u>(2)</u>	1/1000. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

ターゲット情報

機能

Core component of nucleosome. Nucleosomes wrap and compact DNA into chromatin, limiting DNA accessibility to the cellular machineries which require DNA as a template. Histones thereby play a central role in transcription regulation, DNA repair, DNA replication and chromosomal stability. DNA accessibility is regulated via a complex set of post-translational modifications of histones, also called histone code, and nucleosome remodeling.

配列類似性

翻訳後修飾

Belongs to the histone H2A family.

The chromatin-associated form is phosphorylated on Thr-121 during mitosis.

Deiminated on Arg-4 in granulocytes upon calcium entry.

Monoubiquitination of Lys-120 by RING1 and RNF2/RING2 complex gives a specific tag for epigenetic transcriptional repression and participates in X chromosome inactivation of female mammals. It is involved in the initiation of both imprinted and random X inactivation. Ubiquitinated H2A is enriched in inactive X chromosome chromatin. Ubiquitination of H2A functions downstream of methylation of 'Lys-27' of histone H3. Monoubiquitination of Lys-120 by RNF2/RING2 can also be induced by ultraviolet and may be involved in DNA repair. Following DNA double-strand breaks (DSBs), it is ubiquitinated through 'Lys-63' linkage of ubiquitin moieties by the E2 ligase UBE2N and the E3 ligases RNF8 and RNF168, leading to the recruitment of repair proteins to sites of DNA damage. Monoubiquitination and ionizing radiation-induced 'Lys-63'-linked ubiquitination are distinct events.

Phosphorylation on Ser-2 is enhanced during mitosis. Phosphorylation on Ser-2 by RPS6KA5/MSK1 directly represses transcription. Acetylation of H3 inhibits Ser-2 phosphorylation

by RPS6KA5/MSK1.


Symmetric dimethylation on Arg-4 by the PRDM1/PRMT5 complex may play a crucial role in the

germ-cell lineage.

細胞内局在

Nucleus. Chromosome.

画像

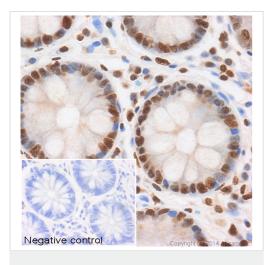

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Histone H2A (acetyl K5) antibody - ChIP Grade (ab1764)

IHC image of ab1764 staining Histone H2A (acetyl K5) in human colon formalin fixed paraffin embedded tissue sections, performed on a Leica Bond. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab1764, 1/1000 dilution, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. No primary antibody was used in the negative control (shown on the inset).

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors