abcam

Product datasheet

Anti-NDRG1 antibody [EPR22461-10] - BSA and Azide free (Capture) ab272294

יעלאעבע RabMAb

画像数 2

製品の概要

製品名

製品の詳細

由来種

アプリケーション

種交差性

免疫原

特記事項

Anti-NDRG1 antibody [EPR22461-10] - BSA and Azide free (Capture)

Rabbit monoclonal [EPR22461-10] to NDRG1 - BSA and Azide free (Capture)

Rabbit

適用あり: Sandwich ELISA

交差種: Human

Recombinant full length protein. This information is proprietary to Abcam and/or its suppliers.

ab272294 is a BSA and Azide Free antibody supplied in an unconjugated format and it is suitable for sandwich ELISAs to quantify Human NDRG1. The recommended pair for sandwich ELISA is:

Capture: ab272294, Human NDRG1 Capture Antibody (unconjugated) Detector: ab272295, Human NDRG1 Detector Antibody (unconjugated)

The reference range value is 390.625 - 25000 pg/ml.

Our carrier-free antibodies are typically supplied in a PBS-only formulation, purified and free of BSA, sodium azide and glycerol. The carrier-free buffer and high concentration allow for increased conjugation efficiency.

This conjugation-ready format is designed for use with fluorochromes, metal isotopes, oligonucleotides, and enzymes, which makes them ideal for antibody labelling, functional and cellbased assays, flow-based assays (e.g. mass cytometry) and Multiplex Imaging applications.

Use our conjugation kits for antibody conjugates that are ready-to-use in as little as 20 minutes with <1 minute hands-on-time and 100% antibody recovery: available for fluorescent dyes, HRP, biotin and gold.

The recommended antibody orientation is based on internal optimization for ELISA-based assays. Antibody orientation is assay dependent and needs to be optimized for each assay type. Please note that the range provided for this antibody is only an estimation based on the performance of the product using the recommended antibody pair. Performance of the antibody pair will depend on the specific characteristics of your assay. We guarantee the product works in sandwich ELISA, but we do not guarantee the sensitivity or dynamic range of the antibody in your assay.

This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility
- Improved sensitivity and specificity
- Long-term security of supply
- Animal-free production

For more information see here.

Our RabMAb[®] technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb**[®] **patents**.

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C.

パッファー Constituent: 100% PBS

キャリア・フリー はい

精製度 Protein A purified

ポリ/モノ モノクローナル **ウローン名** EPR22461-10

アイソタイプ lgG

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab272294の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
Sandwich ELISA		Use at an assay dependent concentration. Can be paired for Sandwich ELISA with PE / R-Phycoerythrin Conjugation Kit - Lightning-Link® (ab102918) and Biotinylation Kit / Biotin Conjugation Kit (Fast and Type A) - Lightning-Link® (ab201795).

ターゲット情報

機能 May have a growth inhibitory role.

組織特異性 Ubiquitous; expressed most prominently in placental membranes and prostate, kidney, small

intestine, and ovary tissues. Reduced expression in adenocarcinomas compared to normal tissues. In colon, prostate and placental membranes, the cells that border the lumen show the

highest expression.

関連疾患 Defects in NDRG1 are the cause of Charcot-Marie-Tooth disease type 4D (CMT4D)


[MIM:601455]; also known as hereditary motor and sensory neuropathy Lom type (HMSNL). CMT4D is a recessive form of Charcot-Marie-Tooth disease, the most common inherited disorder of the peripheral nervous system. Charcot-Marie-Tooth disease is classified in two main groups on the basis of electrophysiologic properties and histopathology: primary peripheral demyelinating neuropathy and primary peripheral axonal neuropathy. Demyelinating CMT

neuropathies are characterized by severely reduced nerve conduction velocities (less than 38 m/sec), segmental demyelination and remyelination with onion bulb formations on nerve biopsy, slowly progressive distal muscle atrophy and weakness, absent deep tendon reflexes, and hollow feet. By convention, autosomal recessive forms of demyelinating Charcot-Marie-Tooth disease are designated CMT4.

配列類似性 細胞内局在 Belongs to the NDRG family.

Cytoplasm. Nucleus. Cell membrane. Whereas in prostate epithelium and placental chorion it is located in both the cytoplasm and the nucleus, nuclear staining is not observed in colon epithelium cells. Instead its localization changes from the cytoplasm to the plasma membrane during differentiation of colon carcinoma cell lines in vitro.

画像

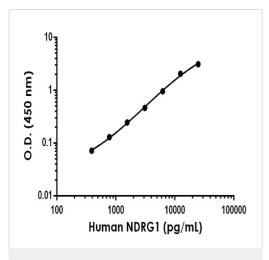

Sandwich ELISA - Anti-NDRG1 antibody [EPR22461-10] - BSA and Azide free (Capture) (ab272294) Representative standard curve from corresponding SimpleStep ELISA® Kit (ab270214).



Anti-NDRG1 antibody [EPR22461-10] - BSA and Azide free (Capture) (ab272294)

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- · Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors