abcam

Product datasheet

Anti-TIMP2 antibody [3A4] ab1828

★★★★★ 7 Abreviews 54 References 画像数 1

製品の概要

製品名 Anti-TIMP2 antibody [3A4]

製品の詳細 Mouse monoclonal [3A4] to TIMP2

由来種 Mouse

 アプリケーション
 適用あり: IHC-P

 種交差性
 交差種: Human

交差が予測される動物種: Mouse, Rat, Rabbit, Guinea pig, Cow

A

免疫原 Synthetic peptide corresponding to Human TIMP2 (N terminal).

ポジティブ・コントロール IHC-P: Normal human placenta tissue.

特記事項 This antibody clone is manufactured by Abcam. If you require a custom buffer formulation or

conjugation for your experiments, please contact orders@abcam.com.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

製品の特性

製品の状態 Liquid

保存方法 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long

term. Avoid freeze / thaw cycle.

ארע"ד Preservative: 0.02% Sodium azide

Constituent: PBS

精製度 Protein G purified

ポリ/モノ モノクローナル

クローン名 3A4 **アイソタイプ** IgG2a

1

アプリケーション

The Abpromise guarantee <u>Abpromise保証は、</u>次のテスト済みアプリケーションにおけるab1828の使用に適用されます アプリケーションノートには、推奨の開始希釈率がありますが、適切な希釈率につきましてはご検討ください。

アプリケーション	Abreviews	特記事項
IHC-P	****(1)	Use a concentration of 1 µg/ml. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

ターゲット情報

機能

Complexes with metalloproteinases (such as collagenases) and irreversibly inactivates them by binding to their catalytic zinc cofactor. Known to act on MMP-1, MMP-2, MMP-3, MMP-7, MMP-8,

MMP-9, MMP-10, MMP-13, MMP-14, MMP-15, MMP-16 and MMP-19.


配列類似性 Belongs to the protease inhibitor l35 (TIMP) family.

Contains 1 NTR domain.

翻訳後修飾 The activity of TIMP2 is dependent on the presence of disulfide bonds.

細胞内局在 Secreted.

画像

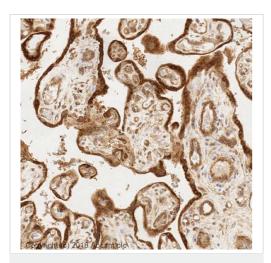

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-TIMP2 antibody [3A4] (ab1828)

IHC image of TIMP2 staining in a section of formalin-fixed paraffinembedded normal human placenta* performed on a Leica BONDTM system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20mins. The section was then incubated with ab1828, 1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.co.jp/abpromise or contact our technical team.

Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors