

# Native Human Fibronectin protein ab80021

画像数 1

### 製品の詳細

製品名	Native Human Fibronectin protein
精製度	> 95 % SDS-PAGE. Purified from human plasma
発現系	Native
アクセッション番号	<b><u>P02751</u></b>
タンパク質長	Full length protein
Animal free	No
由来	Native
生物種	Human
配列の追加情報	Molecular weight: the MW of the hetrodimer is 550 kDa

### 特性

Our **Abpromise guarantee** covers the use of **ab80021** in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

アプリケーション	SDS-PAGE
製品の状態	Lyophilized


### 前処理および保存

保存方法および安定性	Shipped at 4°C. After reconstitution store at -20°C. Avoid freeze / thaw cycles. pH: 7.40 Constituents: 0.476% HEPES, 0.87% Sodium chloride
再構成	Reconstitute with 1.0 milliliters of water to achieve 2.0 mg/mL. Human Fibronectin is known for its solubility issues. The best way to get this material into solution, after reconstitution, is to incubate in a 37°C water bath for at least 30 minutes. Typically this is sufficient, however a small amount of un-dissolved material may remain. If this is the case, a few gentle aspirations with a pipette is recommended.

### 関連情報

<b>機能</b>	<p>Fibronectins bind cell surfaces and various compounds including collagen, fibrin, heparin, DNA, and actin. Fibronectins are involved in cell adhesion, cell motility, opsonization, wound healing, and maintenance of cell shape. Involved in osteoblast compaction through the fibronectin fibrillogenesis cell-mediated matrix assembly process, essential for osteoblast mineralization. Participates in the regulation of type I collagen deposition by osteoblasts.</p> <p>Anastellin binds fibronectin and induces fibril formation. This fibronectin polymer, named superfibronectin, exhibits enhanced adhesive properties. Both anastellin and superfibronectin inhibit tumor growth, angiogenesis and metastasis. Anastellin activates p38 MAPK and inhibits lysophospholipid signaling.</p>
<b>組織特異性</b>	Plasma FN (soluble dimeric form) is secreted by hepatocytes. Cellular FN (dimeric or cross-linked multimeric forms), made by fibroblasts, epithelial and other cell types, is deposited as fibrils in the extracellular matrix. Ugl-Y1, Ugl-Y2 and Ugl-Y3 are found in urine.
<b>関連疾患</b>	Glomerulopathy with fibronectin deposits 2
<b>配列類似性</b>	<p>Contains 12 fibronectin type-I domains.</p> <p>Contains 2 fibronectin type-II domains.</p> <p>Contains 16 fibronectin type-III domains.</p>
<b>発生段階</b>	Ugl-Y1, Ugl-Y2 and Ugl-Y3 are present in the urine from 0 to 17 years of age.
<b>翻訳後修飾</b>	<p>Sulfated.</p> <p>It is not known whether both or only one of Thr-2064 and Thr-2065 are/is glycosylated.</p> <p>Forms covalent cross-links mediated by a transglutaminase, such as F13A or TGM2, between a glutamine and the epsilon-amino group of a lysine residue, forming homopolymers and heteropolymers (e.g. fibrinogen-fibronectin, collagen-fibronectin heteropolymers).</p> <p>Phosphorylated by FAM20C in the extracellular medium.</p> <p>Proteolytic processing produces the C-terminal NC1 peptide, anastellin.</p>
<b>細胞内局在</b>	Secreted, extracellular space, extracellular matrix.

## 画像

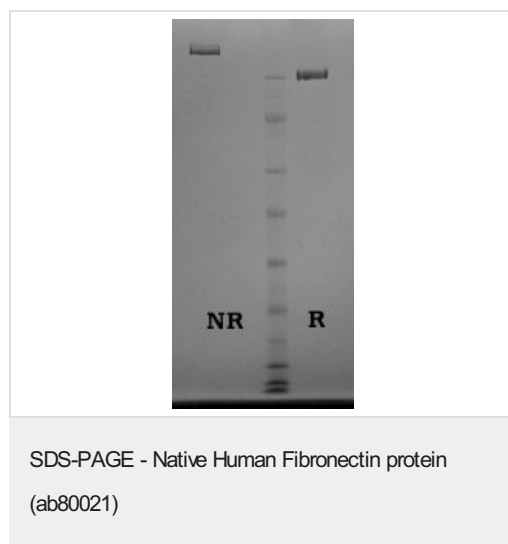

**Gel:** Novex 4-12% Bis-Tris

**Load Load:** Human Fibronectin, 1 µg per lane in Non-reducing (NR) and reducing (R) conditions.

**Buffer:** MOPS

**MW Ladder:** Myosin (191 kDa), Phosphorylase B (97 kDa), BSA (64 kDa), Glutamic Dehydrogenase (51 kDa), Alcohol Dehydrogenase (39 kDa), Carbonic Anhydrase (28 kDa), Myoglobin Red (19 kDa), Lysozyme (14 kDa)

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

### **Our Abpromise to you: Quality guaranteed and expert technical support**

---

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
  
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.co.jp/abpromise> or contact our technical team.

### **Terms and conditions**

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors